meriden

Traveller & Home Irrigator

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www.facebook.com/meridenglobal
www.instagram.com/meridenglobal

Designed by Meriden in Poland
Made in P.R.C.

Model: MD 1522

meriden

Traveller & Home Irrigator


Before using this product, please read carefully this user's manual and keep it properly for later reference.

User manual

MD 1522

Basic introduction

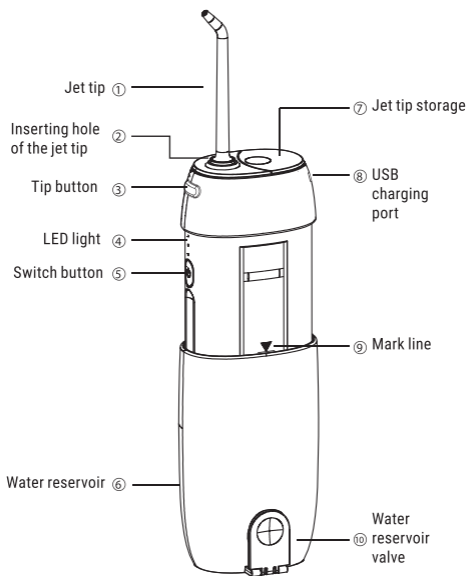Classic jet tip: Removes dental plaque and deeply cleans the teeth.


Working modes: Normal Soft Pulse
Rated voltage: 5V DC
Rated power: 8W
Rated water pressure: 30-90PSI
Pulse frequency: 1600 strokes/min.
Charging time: 3 Hours
Water tank capacity: 120 ML
Battery Spec.: 1500mAh 3.7V Lithium battery
Executive standard: GB 4706.1-2005 GB 4706.59-2008
Q/RMZN 5-2017

User instruction


Charging

- a. Switch the device off.
- b. Push the water reservoir back to its original position.
- c. Connect the Micro USB charger to the water flosser charging port.
- d. Connect the adapter plug to power supply and make sure that the indicator light is on.
- e. When the device is being charged, the indicator is flashing white and stops flashing after being fully charged. The device can be used for 60 minutes after a single full-time charging for 3 hours.

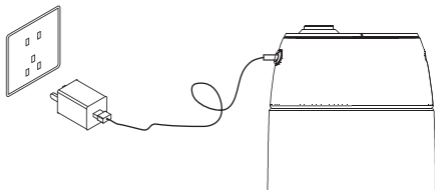Installing and removing the jet tip

- Take out classic jet tip from the tip storage.
- Insert the classic jet tip into hole until the tip is held firmly.
- A 'click' sound will mean that the jet tip is installed.
- When removing the jet tip, press the tip button first, then take the jet tip out and put it into the tip storage.

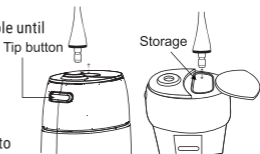
Warning: Please do not press the tip button when the device is working.

Water filling method

- Open the valve of the water reservoir and pull the water reservoir downward toward the water line.
- Hold the water reservoir horizontally and pour water into the water reservoir.
- Close the valve of the water reservoir.

Note: Before pulling the water reservoir downward, open the water reservoir valve, otherwise the air pressure will make pulling the water reservoir difficult.

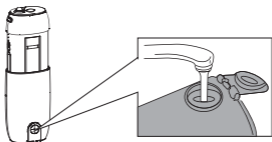
Running water is recommended for cleaning purposes. If you feel the water is too cold, you can use warm water with temperature below 40°C as well as salty water with concentration below 50% or mouthwash.

Modes and the functional switch

- Press the switch button, the indicator light is on.
- Change between modes – 'Normal', 'Soft', 'Pulse' by pressing the switch button within 5 seconds.

The 'Soft' mode is recommended for the first use of the device. The product is equipped with smart memory chip which memorizes your habits. When switched on, the device will automatically work on the mode you used the last time.

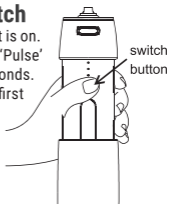Operating instructions:

- Put the jet tip into your mouth and press the switch button to select the desired mode. Aim the jet tip vertically at the teeth, crevices or gums to ensure that the area is rinsed properly.
- Pay close attention to areas between the teeth and around the tooth socket.
- While cleaning your teeth, open your mouth slightly to ensure that water flows out.
- Keep the device vertically during use to prevent water flow from being stopped by detached water pipe inside the water reservoir.

Cleaning and Maintenance

Cleaning

- Switch off the power, open water reservoir cover and pour the water out.
- Turn on the switch to drain the water from reservoir.
- The product has IPX7 waterproof design and can be directly cleaned with running water.
- After thorough cleaning, drain the water and leave the reservoir to dry.

Maintenance

Tip

- Wash with water and dry with a cloth.
- It is recommended to replace the tip every 6 months. If the tip is deformed, please change it in time, otherwise deformed tip will affect the use and may cause the device to break down.

Body

- Use a soft cloth lightly moistened with small amount of detergent to clean the outside of the product.
- The water reservoir unit is removable for easy cleaning.
- It is recommended to clean water reservoir once a week.

Note:

- Do not wash the device with hot water above 50°C.
- Do not use diluents, alcohol or other stimulants as they will discolor or damage the device.

Attention

To prevent burning, electrocution, fire or injury to people, take the following safety precautions:

- Do not plug this device into a voltage system different from the voltage system specified on the device or charger.
- Do not touch the charger with wet hands.
- Do not use the product if there is any damage to the cord, plug or the product itself or the product is not working properly.
- Do not direct water under the tongue, into an ear, nose or other delicate areas, as the pressure may injure them seriously.
- Use this product only as indicated in these instructions or as recommended by your dentist.
- Fill the reservoir with water or dentist-recommended solutions only.
- Be careful not to drop the product. Never push any objects through the openings in this product. Power cord cannot be replaced- in the case of the power cord damage, the product must not be used.
- Do not use this product outdoors or operate it in the area where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use iodine, bleach or water insoluble concentrated essential oils in this product. The use of such substances may lead to impairment of the product performance and shortening of the product life.
- Remove any oral jewellery before using this product.
- Do not use this product if you have any open wounds on the tongue or in the mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Children must be supervised to ensure that they do not play with the device.
- Storage place should be dry and clean, as well as kept away from direct sunlight and heat.
- If the device needs to be stored for a period of time longer than 3 months, it must be fully charged before storing.
- It is recommended not to operate the device on an empty water reservoir in order to avoid the product damage.

Environmental Protection

If you want to discard this product, please, do not place it together with domestic waste, but send it to the officially designated recycling center and dispose of the battery at the officially designated recycling center. Doing so can support environmental protection.

Dealing with the battery

The built-in rechargeable battery contains substances which may pollute the environment.

Please, make sure to remove the built-in battery only when you intend to dispose of the device. Then, take the battery and the device to the official recycling center.

Do not remove the rechargeable battery of the device if you do not intend to dispose of it. Make sure the battery has been exhausted before removing it.

Removing the battery

1. Switch off the device.
2. Open the device top cover.
3. Remove the battery.

Troubleshooting

Problem	Reasons	Solutions
The device does not work	The product is new - it has just been purchased or it has been idle for more than 3 months	Connect to the charging base, charge the device for 24 hours.
There is leakage between the brush head and the handle knob	The brush head is not completely inserted into the handle knob	Take off the brush head and replace it with a new one.
Water pressure is low	Dirt has stuck in the tip. The top of the tip has been deformed. Dirt has got into the pump.	Replace the tip. Contact the customer service.
Water does not come out	Dirt has stuck in the filter. Water reservoir is not filled with water. The water irrigator has been tilted too much when used. The water pipe has not been installed. The water pump is full of dirt.	Clean the filter. Fill the reservoir with water. Use the device in the vertical position. Install the pipe. Contact the customer service

Meriden Sp. z o.o.
ul. Piastów 96
42-290 Blachownia, Poland
NIP: 9492217625
meriden@meriden.pl
www.meridenglobal.com

www.facebook.com/meridenglobal
www.instagram.com/meridenglobal

Designed by Meriden in Poland
Made in P.R.C.

Model: MD 1522