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Welcome

Thank you for purchasing our product. We recommend you to use it with the supporting app to learn more about maintaining a healthy smile.

Important Safety Tips

It is recommended that this toothbrush be used only by a person over 12 years old.

It is normal to feel slightly uncomfortable when the toothbrush is used for the first time.

The toothbrush can only be used with Oclean special heads. We will not bear any responsibility due to the use of any third-party accessory.

This toothbrush is intended for personal brushing only. It is not recommended that the toothbrush be shared with others.

Avoid exposing the toothbrush to direct sunlight. Please power off the charging holder when it is not in use and before cleaning it.

Not recommended for use

Individuals who have recently undergone oral surgery. It is recommended that consumers with oral disease, such as dental cavities, periodontal (gum) disease, oral cancer, oral infectious diseases, consult with a dentist before using this toothbrush.

Individuals with mental and/or physical disabilities should only be used under supervision.

People who have mouth diseases, loose false teeth or gums with hemorrhage.

People with sensitive and susceptible oral cavity, they should consult a professional dentist before using this product.

The toothbrush is not suitable for people with reduced mental and physical capabilities. Not suitable for infant and children under the age of 12. Children who use this product should be strictly supervised by an adult.

Instructions

First Time Use

The toothbrush can be turned on by pressing the ON/OFF button for the first time on a brand-new toothbrush. The bluetooth function will be activated and remain on. After installing the Oclean App, please complete the prompts to complete the connection and settings for the toothbrush.

Install App & Toothbrush Settings

You can use the app to customize your brushing programs - automating tooth cleaning, tooth whitening and gum massaging in one brushing session. There is no need to open the Oclean App each time you brush. Simply open the app every few days to sync the brushing data.

Direct Use without App Settings


You can also use the toothbrush without connecting to the app. It has a touchscreen display and an ON/OFF button. It offers touch control via the touch-sensitive screen.

Awake and Automatic Sleep


The toothbrush can be turned on by pressing the ON/OFF button. When the toothbrush is not in use for 15s, it will automatically enter the sleep mode.

Default Brushing Plan & Direct Brushing


Press the ON/OFF button again after the toothbrush is awakened, it will start working with the default "Standard Cleaning" and lasts for two minutes.

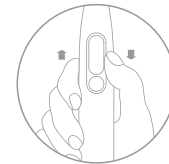Slide up/down

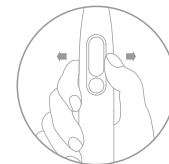Slide right/left

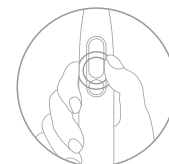Long press on screen or short press on power button to confirm the settings

Brushing Plan & App Settings

You can either choose or customize a brushing plan through the supporting Oclean App. You can also upgrade the firmware to continuously update its performance.

Low Battery Prompt

When battery percentage is lower than 20%, you'll see low battery prompt appears on your toothbrush display frequently, you will need to charge the toothbrush. When the battery percentage is lower than 1%, you will be prompted by motor vibration, and the screen will not be on.

How to Charge

Remove the charging base from the wall mount and connect it to the mobile phone charger. When the bottom of the toothbrush is placed gently above the charging base, it will automatically position itself. The display will show the battery percentage when the toothbrush is charging. A full charge usually takes about two hours. (Recommended adapter: 5V/1A - 5V/2A).

Important Charging Tips

The toothbrush runs with a fast-charge mechanism. It is recommended that the toothbrush be placed on the charging base only when it is necessary. Please do not put the toothbrush on the charging base for a long time, as it may cause the built-in battery to decay prematurely so that the user experience is affected.

- 1. No disassembly, impact, extrusion or throwing into fire
- 2. Please stop using the toothbrush if the battery has expanded
- 3. Please don't put the battery in high temperature environment
- 4. Do not use the battery if water entered the device.

• Specifications

Product Name	Oclean Smart Sonic Electric Toothbrush
Model	Oclean X
Version	Pro
Max Torsion	220gf.cm
Frequency	2402-2480MHz
Max Output	0dBm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V = 1A
Battery Life*	20 days (Level 1), 40 days (Level 32)
Battery Capacity	800 mAh
Charging Time	2 hours
Bluetooth	4.2 above
Standard	GB 39669-2020

* Calculated by 2 minutes/time, twice/day.

EU REP

EU Authorized Representative
 Agency name: Mbrush UAB
 Representative: Algirdas Tevelis
 European Address: Totorių g.3, Kaunas, Lithuania
 Post Code: LT44237
 Email: service@mobalt.eu
 Phone: +37067052772

Facebook Page

www.facebook.com/ocleanofficial

Follow the Oclean Official Facebook Page

Download App

www.oclean.com/download

Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
- × The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- * Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

Name of parts	Toxic and harmful substance or element						Environment Protection Term
	Pb	Hg	Cd	Cr6+	PBB	PBDE	
Vibrating motor	○	○	○	○	○	○	10 Years
Accessories	○	○	○	○	○	○	5 Years
Battery	○	○	○	○	○	○	

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Warranty Card

Delivery Date	Employed equipment	Acceptance date	Service man

Please ask for a maintenance report from the maintenance center and properly keep it with the warranty certificate. Shenzhen Yunding Information Technology Co., Ltd., in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers guarantee for customers who purchase Oclean Smart Sonic Electric Toothbrush. In case of inconsistency between the service articles and the provisions of guarantee, the provisions of guarantee and the contents in the guarantee certificate shall govern.

Guarantee Certificate of Oclean Smart Sonic Electric Toothbrush

User Information

Name: _____ Tel: _____ E-Mail: _____

Address: _____ Zip Code: _____

Sales information

Name of Retailer: _____ Tel: _____

Address: _____ Zip Code: _____

Invoice No: _____ Sales Date: _____

Type: _____

Production Place: _____

SN

FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.
- The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

FCC ID: 2AN5D-OCLEANX

IC WARNING

This device contains licence-exempt transmitter (s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:
 (1) This device may not cause interference.
 (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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Made in China

If the product you have purchased is found to have a fault or defect, please contact us.
 North America
 Call us at - (866)OCLEAN9 or (866) 625-3269
 Email: support.us@oclean.com

Other Region
 Email: service@oclean.com

Manufacturer: Shenzhen Yunding Information Technology Co., Ltd.

Manufacturer's address: #28G, Building 3, Dachong Business Center (west of Tonggu Road and north of Shennan Avenue), Yuehai Subdistrict, Nanshan District, Shenzhen, China

Hereby, Shenzhen Yunding Information Technology Co., Ltd. declares that the radio equipment type Y3 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: <http://www.oclean.com/en/service.aspx>

For more regulatory information and product certification/compliance logos related to Oclean X Pro, please go to Settings->About->Regulatory on your toothbrush.

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Producent: Shenzhen Yunding Information Technology Co., Ltd.