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Oclean Sonic Electric Toothbrush User Manual


Thanks for using Oclean Sonic Electric Toothbrush.

Important Safety Notice

- When using the product for the first time, due to the high-speed swinging of brush, if it feels stimulating or numb, please do not worry, it is normal and the feeling will weaken and disappear gradually in around 7-14 days.
- If not using for long time (the product battery has no power due to the standby power consumption), when charging for the toothbrush again, the charging indicator doesn't show, please wait for 6-10 minutes patiently, it can recover normal charging.
- If the main part of toothbrush becomes slightly hot when charging, it is normal. If you feels hot with your hand, please immediately stop charging and submit the toothbrush to Oclean's After-sales Services Center for testing.
- The product can only use the special brush head for Oclean Sonic electric toothbrush. Oclean will not be held liable for all the consequences caused by the use of other unofficial or non-licensed brush head products.
- The product is only used for cleaning teeth and must not be used for other purposes.
- The product is only intended as a personal cleaning product. It is not recommended to use it by multiple persons simultaneously.
- The product is only used indoors, do not place the toothbrush under the blazing sun for quite a long time. It shall be placed in a dry and well-ventilated place with normal temperature.

Attention should be paid for those with the following conditions of oralsensitivities

- Those who have recently undergone periodontal surgery or are undergoing treatments.
- Individuals with other symptoms in the mouth including those with unstable dentures or have persistent heavy bleeding from the gums.
- Other individuals with sensitive and fragile oral issues. This product is not recommended for use with the above group or consult a professional dentist before use.
- This product is not suitable for people with disabilities, mental illness or lack of relevant experience and knowledge (including children (0-12 years old)). Unless used under the supervision of their guardians.

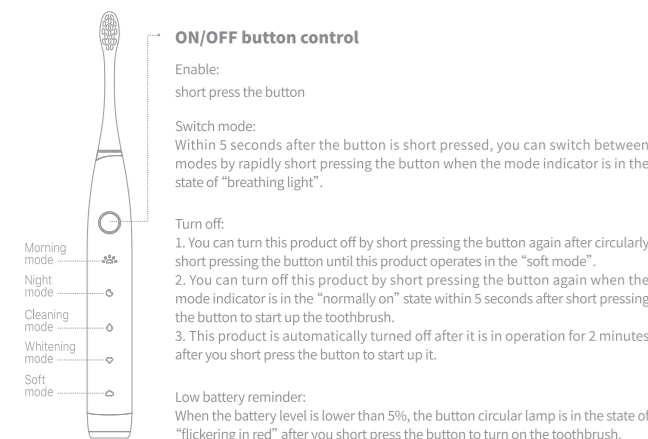If you need more detailed product usage information, please follow the Oclean Official Account.

Instruction for use

Default Brushing program

Five toothbrushing modes are built in this product: Morning mode, night mode, Cleaning mode, whitening mode and soft mode.

- Morning mode: compared with "night mode", the pressure is low. As a cleaning and care tool, it can effectively remove dental plaque caused during the sleep.
- Night mode: deeply remove food residues and dental plaque on the teeth with high pressure.
- Cleaning mode: meet requirements for routine and ordinary cleaning. Compared with "night mode", the intensity level would be higher.
- Whitening mode: effectively clean stains on those who have habits of smoking or drinking.
- Soft mode: protect the health of gingiva while fulfilling the fundamental function of removing dental plaque in those who have oral allergy.

Type-C direct charging

Open the bottom cover of the electric toothbrush, connect Type-C interface into the port at the bottom of the toothbrush, connect one end of the USB interface to the power supply, the button circular light would flash in white and this indicates that the toothbrush is being charged.

Charging process

Charging status: the button circular light is on white and flash

Use for the first time

The user may charge the toothbrush depending on his or her need for battery life to ensure normal use of this product.

Application process

When the battery level is lower than 5%, the button circular light is in the state of "flickering in red" after you short press the button to turn on the toothbrush.

Precautions

Battery can not be replaced by users themselves, it can only be replaced by the personnel in the customer center authorized by Oclean, non-specialized personnel can not dismantle the product by themselves, or it will be regarded as that the user give up the warranty automatically.

Please charge the Oclean toothbrush with a 5V/1A or 5V/2A mobile phone charger.

Such mobile phone charger includes iPhone charger, tablet charger, etc. If the charging current is lower than 1A, this may cause the charging time to be extended; if you charge the toothbrush with a mobile phone charger with a voltage of higher than 5V, this may cause the charger circuit to be damaged.

Recommended tooth brushing method

Five toothbrushing modes are built in this product: Morning mode, night mode, Cleaning mode, whitening mode and soft mode.

Move the toothbrush slowly at 45° between the brush and gingiva line. There is no need to manually brush up and down during this process.

Hold the toothbrush handle gently, move slowly in a circular motion, and please do not apply excessive pressure.


Please gargle with clear water after brushing.

Washing and maintenance

Cleaning the toothbrush body and the brush head

After rinsing with tap water, put on the cover and wait for the brush head to dry naturally, do not pull the bristles of the brush head.

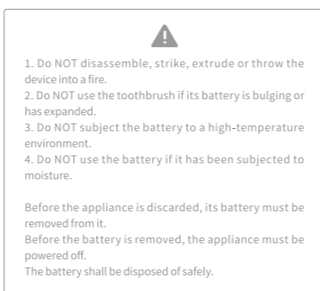Product specifications

Product Name	Oclean Sonic Electric Toothbrush
Model	F5002
IP Rating	IPX7
Power of Subject	5W
Charger	Type-C Direct Charging
Input	5V/1A
Service Time*	≥ 180 Days (**)
Battery Capacity	2500 mAh(**)
Charging time	≤ 5 hours(**)
Implementation standard	GB 4706.1-2005 GB 4706.59-2008

* Validated by 2 brushing processes per day and 2 minutes of each process.

** The battery capacity theoretically has a capacity tolerance of ±50mAh.

*** It is suggested to use a standard 5V/1A adapter.

Please request a maintenance report from the maintenance center and keep it together with the warranty certificate.

Shenzhen Yunding Information Technology Co., Ltd., in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers a guarantee for consumers who purchase the Oclean Smart Sonic Electric Toothbrush. In the case of any inconsistency between the service articles and the provisions of guarantee, the provisions of the guarantee and the contents of the guarantee certificate shall govern.

If the product you have purchased is found to have a fault or defect, please contact our service line at 4008-365-909.

Manufacturer: Shenzhen Yunding Information Technology Co., Ltd.

Manufacturer's address: #28G, Building 3, Dachong Business Center (west of Tonggu Road and north of Shennan Avenue), Yuehai Subdistrict, Nanshan District, Shenzhen, China

Website: www.oclean.com

Remarks
This indicates that the content of toxic and hazardous substances in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011.
This indicates that the content of toxic and hazardous substances in at least one of the homogeneous materials of this component is more than the limited quantity as specified in GB/T 26572-2011.
Except for very few materials, due to technical development reasons, most of the materials used are environmentally friendly. Please use this product with confidence during the environmental protection period.

Name of parts	Toxic and harmful substance or element				Environment Protection	10 years
	DEHP/BBP/DBP/DIBP	PBDEs	PBBs	Cr6+		
Vibrating motor	○	○	○	○	○	○
Accessories	○	○	○	○	○	○
Battery	○	○	○	○	○	○

Name of parts	Toxic and harmful substance or element				Environment Protection	10 years
	DEHP/BBP/DBP/DIBP	PBDEs	PBBs	Cr6+		
PCBA	○	○	○	○	○	○
Plastic parts	○	○	○	○	○	○
Metal hardware	○	○	○	○	○	○

Guarantee Certificate of Oclean Sonic Electric Toothbrush

Product Information

Name: _____ Phone: _____ E-mail: _____ Invoice no.: _____ Product type: _____
Address: _____ ZIP code: _____ Retailer: _____ Address: _____ Sales date: _____ Manufacturing site: _____

Warranty Card of Oclean Sonic Electric Toothbrush

Delivery date	Equipment employed	Service date

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