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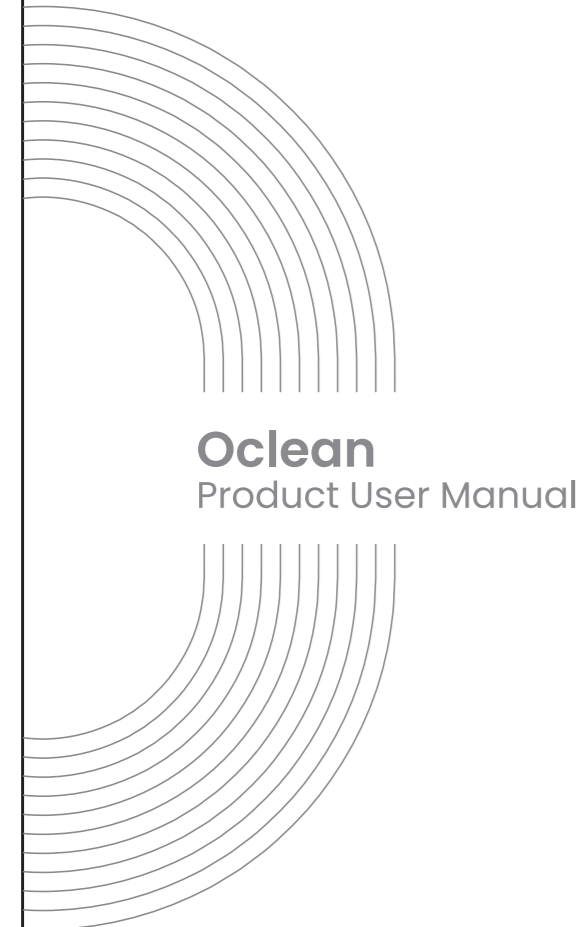
### Oclean Sonic Electric Toothbrush User Manual


Thanks for using Oclean Sonic Electric Toothbrush.

#### Important Safety Notice

- When using the product for the first time, due to the high-speed swinging of brush, if it feels stimulating or numb, please do not worry, it is normal and the feeling will weaken and disappear gradually in around 7-14 days.
- If not using for long time (the product battery has no power due to the standby power consumption), when charging for the toothbrush again, the charging indicator doesn't show, please wait for 6-10 minutes patiently, it can recover normal charging.
- If the main part of toothbrush becomes slightly hot when charging, it is normal. If you feels hot with your hand, please immediately stop charging and submit the toothbrush to Oclean's After-sales Services Center for testing.
- The product can only use the special brush head for Oclean Sonic electric toothbrush. Oclean will not be held liable for all the consequences caused by the use of other unofficial or non-licensed brush head products.
- The product is only used for cleaning teeth and must not be used for other purposes.
- The product is only intended as a personal cleaning product. It is not recommended to use it by multiple persons simultaneously.
- The product is only used indoors, do not place the toothbrush under the blazing sun for quite a long time. It shall be placed in a dry and well-ventilated place with normal temperature.

#### Attention should be paid for those with the following conditions of oral sensitivities

- Those who have recently undergone periodontal surgery or are undergoing treatments.
- Individuals with other symptoms in the mouth including those with unstable dentures or have persistent heavy bleeding from the gums.
- Other individuals with sensitive and fragile oral issues. This product is not recommended for use with the above group or consult a professional dentist before use.
- This product is not suitable for people with disabilities, mental illness or lack of relevant experience and knowledge (including children (0-12 years old)). Unless used under the supervision of their guardians.

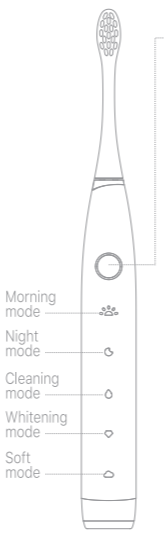**ON/OFF button control**  
 Enable: short press the button  
 Switch mode: Within 5 seconds after the button is short pressed, you can switch between modes by rapidly short pressing the button when the mode indicator is in the state of "breathing light".  
 Turn off: 1. You can turn this product off by short pressing the button again after circularly short pressing the button until this product operates in the "soft mode".  
 2. You can turn off this product by short pressing the button again when the mode indicator is in the "normally on" state within 5 seconds after short pressing the button to start up the toothbrush.  
 3. This product is automatically turned off after it is in operation for 2 minutes after you short press the button to start up it.  
 Low battery reminder: When the battery level is lower than 10%, the button circular lamp is in the state of "flickering in red" after you short press the button to turn on the toothbrush.

#### Instruction for use

- Default Brushing program**  
 Five toothbrushing modes are built in this product: Morning mode, night mode, cleaning mode, whitening mode and soft mode.
- Morning mode:** compared with "night mode", the pressure is low. As a cleaning and care tool, it can effectively remove dental plaque caused during the sleep.
  - Night mode:** deeply remove food residues and dental plaque on the teeth with high pressure.
  - Cleaning mode:** meet requirements for routine and ordinary cleaning. Compared with "night mode", the pressure is relatively tender.
  - Whitening mode:** effectively clean bite marks in those who have habits of smoking or frequently drink coffee, tea, etc.
  - Soft mode:** protect the health of gingiva while fulfilling the fundamental function of removing dental plaque in those who have oral allergy.

#### Type-C direct connection charging

Open the bottom cover of the electric toothbrush, connect Type-C interface into the port at the bottom of the toothbrush, connect one end of the USB interface to the power supply, the button circular lamp flickers in white and this indicates that the toothbrush is being charged.

#### Charging process

Charging status: the button circular lamp is on in white and flickers.

#### Use for the first time

The user may charge the toothbrush depending on his or her need for battery life to ensure normal use of this product.

#### Application process

When the battery level is lower than 10%, the button circular lamp is in the state of "flickering in red" after you short press the button to turn on the toothbrush.

#### Precautions

Battery can not be replaced by users themselves, it can only be replaced by the personnel in the customer center authorized by Oclean, non-specialized personnel can not dismantle the product by themselves, or it will be regarded as that the user give up the warranty automatically.


Please charge the Oclean toothbrush with a 5V/1A or 5V/2A mobile phone charger.

Such mobile phone charger includes iPhone charger, tablet charger, etc. If the charging current is lower than 1A, this may cause the charging time to be extended; if you charge the toothbrush with a mobile phone charger with a voltage of higher than 5V, this may cause the charger circuit to be damaged.


**1. Do NOT** disassemble, strike, extrude or throw the device into a fire.  
**2. Do NOT** use the toothbrush if its battery is bulging or has expanded.  
**3. Do NOT** subject the battery to a high-temperature environment.  
**4. Do NOT** use the battery if it has been subjected to moisture.

Before the appliance is discarded, its battery must be removed from it. Before the battery is removed, the appliance must be powered off. The battery shall be disposed of safely.

#### Recommended tooth brushing method

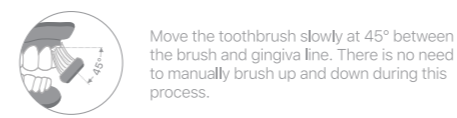Move the toothbrush slowly at 45° between the brush and gingiva line. There is no need to manually brush up and down during this process.

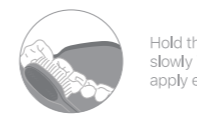
Hold the toothbrush handle gently, move slowly in a circular motion, and please do not apply excessive pressure.

Please gargle with clear water after brushing.

#### Cleaning and Maintenance

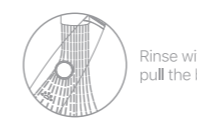Rinse with tap water and air dry, and do not pull the bristles of the brush head

Product specifications

Product Name	Oclean Sonic Electric Toothbrush
Model	F5002
IP Rating	IPX7
Power of Subject	1.85W
Charger	Type-C direct Recharge
Input	5V=1A
Service Time*	≤ 180 Days (*)
Battery Capacity	2500 mAh(**)
Charging time	≤ 5 hours(**)
Implementation standard	Q/JSZYD-1-2021

\* Validated by 2 brushing processes per day and 2 minutes of each process.(Morning mode)  
 \*\* The battery capacity theoretically has a capacity tolerance of ±50mAh.  
 \*\*\*It is suggested to adopt the USB charger of 5V=1A in the phone industry.

#### FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:  
 -- Reorient or relocate the receiving antenna.  
 -- Increase the separation between the equipment and receiver.  
 -- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.  
 -- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

#### IC WARNING

This device contains licence-exempt transmitter (s)/receiver (s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:  
 (1) This device may not cause interference.  
 (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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CAN ICES-3(B)/NMB-3(B)

Hereby, Shenzhen Yunding Information Technology Co., Ltd. declares that the radio equipment type F5002 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: <http://www.oclean.com/en/service.aspx>

Please request a maintenance report from the maintenance center and keep it together with the warranty certificate.

Shenzhen Yunding Information Technology Co., Ltd., in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers a guarantee for consumers who purchase the Oclean Sonic Electric Toothbrush. In the case of any inconsistency between the service articles and the provisions of guarantee, the provisions of the guarantee and the contents of the guarantee certificate shall govern.

If the product you have purchased is found to have a fault or defect, please contact our service line at [service@oclean.com](mailto:service@oclean.com)

Manufacturer: Shenzhen Yunding Information Technology Co., Ltd. Manufacturer's address: #28G, Building 3, Dachong Business Center (west of Tonggu Road and north of Shennan Avenue), Yuehai Subdistrict, Nanshan District, Shenzhen, China  
 Production address: Building A, No.6 of Xinmu Road, Xinmu Community, Pinghu, Shenzhen, China  
[www.oclean.com](http://www.oclean.com)

If you need more detailed product usage information, please follow the Oclean Official Account

[www.facebook.com/ocleanglobal](https://www.facebook.com/ocleanglobal)

**Remarks**  
 This indicates that the content of toxic and hazardous substances in all homogeneous materials of this component is less than the limits specified in the following table:  
 Directive 2011/65/EU & (EU) 2015/863. This indicates that the content of toxic and hazardous substances in at least one of the materials used in this component is more than the limited quantity as specified in RoHS 2.0 Directive (EU) 2015/863. This component is exempt from the scope of the RoHS 2.0 Directive (EU) 2015/863 and amendment 2015/146/ EU. Except for very few materials, due to the provisions of the RoHS 2.0 Directive, the materials used are environmentally friendly. Please use this product with confidence during the environmental protection period.

Name of parts	Toxic and harmful substance or element				Environment Protection	10 years
	DEHP/BBP/DBP/BBP	PBDEs	PBBs	Cd		
Vibrating motor	○	○	○	○	○	○
Accessories	○	○	○	○	○	○
Battery	○	○	○	○	○	○

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Guarantee Certificate of Oclean Sonic Electric Toothbrush

Warranty Card of Oclean Sonic Electric Toothbrush

Name: \_\_\_\_\_ Phone: \_\_\_\_\_ E-mail: \_\_\_\_\_ Address: \_\_\_\_\_ Invoice no.: \_\_\_\_\_ Product type: \_\_\_\_\_  
 Address: \_\_\_\_\_ Phone: \_\_\_\_\_ E-mail: \_\_\_\_\_ Address: \_\_\_\_\_ Sales date: \_\_\_\_\_ ZIP code: \_\_\_\_\_

Manufacturer's address: #28G, Building 3, Dachong Business Center (west of Tonggu Road and north of Shennan Avenue), Yuehai Subdistrict, Nanshan District, Shenzhen, China

Production address: Building A, No.6 of Xinmu Road, Xinmu Community, Pinghu, Shenzhen, China

[www.oclean.com](http://www.oclean.com)

Delivery date	Equipment employed	Acceptance date	Service man

Fault description/reason (disposal)

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Produttore: Shenzhen Yunding Information Technology Co., Ltd.  
Indirizzo del produttore: #28G, Building 3, Dachong Business Center (west of Tonggu Road and north of Shennan Avenue), Yuehai Subdistrict, Nanshan District, Shenzhen, China  
Indirizzo di produzione: Building A, No. 6, Xinmu Avenue, Xinmu Community, Pinghu Street, Longgang District, Shenzhen  
Site web: [www.oclean.com](http://www.oclean.com)

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