Panasonic

Operating Instructions (Household use) Sonic Vibration Toothbrush

Model No. EW-DP52

| | | | | | |
|-------------------|------------|------------------|------------|------------------|------------|
| <u>English</u> | <u>3</u> | <u>Português</u> | <u>135</u> | <u>Česky</u> | <u>267</u> |
| <u>Deutsch</u> | <u>25</u> | <u>Svenska</u> | <u>157</u> | <u>Slovensky</u> | <u>289</u> |
| <u>Français</u> | <u>47</u> | <u>Dansk</u> | <u>179</u> | <u>Magyar</u> | <u>311</u> |
| <u>Italiano</u> | <u>69</u> | <u>Norsk</u> | <u>201</u> | <u>Română</u> | <u>333</u> |
| <u>Español</u> | <u>91</u> | <u>Suomi</u> | <u>223</u> | <u>Türkçe</u> | <u>355</u> |
| <u>Nederlands</u> | <u>113</u> | <u>Polski</u> | <u>245</u> | | |

Operating Instructions
(Household use) Sonic Vibration Toothbrush
Model No. EW-DP52

Contents

| | | | |
|----------------------------------|----------|---|-----------|
| Safety precautions..... | 6 | Changing the brush head (Sold separately)..... | 20 |
| Parts identification | 11 | FAQ | 20 |
| Before cleaning your teeth | 12 | Troubleshooting | 21 |
| Features..... | 12 | Battery life | 23 |
| Charging the toothbrush | 13 | Removing the built-in rechargeable battery..... | 23 |
| Basic usage..... | 14 | Specifications | 24 |
| Tips on usage | 16 | | |
| After use | 18 | | |

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

WARNING

- Be sure to switch off the appliance and remove the power plug from a household outlet before performing any maintenance. Failure to do so may result in electric shock.
- Do not wash the AC adaptor and the charging stand. Doing so may result in electric shock.
- The supply cord cannot be replaced. If the cord is damaged, the AC adaptor should be scrapped.
- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.

- Do not use anything other than the supplied AC adaptor and the charging stand for any purpose. Also, do not use any other product with the AC adaptor and the charging stand. (See page 11.)
- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type references of power supply unit and charging stand are marked near the symbol.

Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

DANGER

Denotes a potential hazard that will result in serious injury or death.

WARNING

Denotes a potential hazard that could result in serious injury or death.

CAUTION

Denotes a hazard that could result in minor injury or property damage.

The following symbols are used to classify and describe the type of instructions to be observed.

This symbol is used to alert users to a specific operating procedure that must not be performed.

This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

WARNING

► Power supply

Do not connect or disconnect the adaptor to a household outlet with a wet hand.

-Doing so may cause electric shock or injury.

Do not immerse the AC adaptor and the charging stand in water or wash it with water. Do not place the AC adaptor and the charging stand over or near water filled sink or bathtub. Do not use when the AC adaptor is damaged or when the fitting into the household outlet is loose.

Do not damage the cord or the power plug.

- Do not deface, modify, forcefully bend or pull the cord or the power plug. Also, do not place anything heavy on, twist or pinch the cord.

-Doing so may cause electric shock or fire due to a short circuit.

WARNING

Do not use in a way exceeding the rating of the household outlet or the wiring.

- Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.

Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the AC adaptor.

Fully insert the adaptor and the appliance plug.
- Failure to do so may cause fire or electric shock.

Always unplug the adaptor from a household outlet when cleaning the AC adaptor and the charging stand.

- Failure to do so may cause electric shock or injury.

Regularly clean the power plug and the appliance plug to prevent dust from accumulating.

- Failure to do so may cause fire due to insulation failure caused by humidity.

Disconnect the adaptor and wipe with a dry cloth.

► In case of an abnormality or malfunction

Immediately stop using and remove the adaptor if there is an abnormality or failure.

- Failure to do so may cause fire, electric shock, or injury.

<Abnormality or failure cases>

• **The main unit, adaptor, charging stand or cord is deformed or abnormally hot.**

• **The main unit, adaptor, charging stand or cord smells of burning.**

• **There is abnormal sound during use or charging of the main unit, adaptor, charging stand or cord.**

- Immediately request inspection or repair at an authorized service centre.

► This product

This product has a built-in rechargeable battery.

Do not throw into fire, apply heat, or charge, use, or leave in a high temperature environment.

- Doing so may cause overheating, ignition, or explosion.

Do not modify or repair.

- Doing so may cause fire, electric shock, or injury.

Contact an authorized service centre for repair (battery change, etc.).

Never disassemble except when disposing of the product.

- Doing so may cause fire, electric shock, or injury.

WARNING

► Preventing accidents

-  **Do not store within the reach of children or infants. Do not let them use it.**
- Doing so may cause an accident or injury due to accidental ingestion of removable parts.

CAUTION

► Preventing damage to teeth and gums

- Do not use the brush with toothpaste that contains a high abrasive content, such as toothpaste for removing tar or for whitening. Do not apply the brush head strongly against teeth or gums, or use it continuously against same spot.**
-  **Do not apply the plastic part of the main unit or brush head to teeth.**
Do not let anyone who cannot operate the brush by themselves use it.
People with low sensitivity in their mouth should not use the brush.
- Doing so may cause injury or damage to teeth and gums.

 **People who may have periodontal disease, have teeth being treated, or who are concerned about symptoms in their mouth should consult a dentist before use.**

- Failure to do so may cause injury or damage to teeth and gums.

 **People who develop allergy symptoms such as itchiness, redness, or rash while using the silicone brush should stop using it immediately and consult a physician.**

- Failure to do so may make symptoms worse.

► This product

 **Never use for a purpose other than cleaning your teeth.**

- Doing so may cause an accident or injury.

 **Do not share the brush head with your family or other people.**

- Doing so may result in infection or inflammation.

 **Do not place a metallic object such as a coin or a clip on top of the charging stand.**

- Doing so may cause burns due to excessive heat.

► Note the following precautions

 **Do not drop or hit against other objects.**

- Doing so may cause injury.

CAUTION

Do not wrap the cord around the adaptor when storing.

- Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.

Do not allow metal objects or trash to adhere to the power plug or the appliance plug.

- Doing so may cause electric shock or fire due to a short circuit.

Disconnect the adaptor from the household outlet when not charging.

- Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

Disconnect the adaptor or the appliance plug while holding onto the adaptor or the appliance plug instead of the cord.

- Failure to do so may cause electric shock or injury.

Persons using medical pacemakers, etc. should consult their physician before use.

- Failure to do so may affect pacemakers, etc., due to electromagnetic waves generated by the charging stand while it is powered.

Handling of the removed battery when disposing

DANGER

- **The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.**
 - **Do not charge the battery after it has been removed from the product.**
 - **Do not throw into fire or apply heat.**
 - **Do not hit, disassemble, modify, or puncture with a nail.**
 - **Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.**
 - **Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.**
 - **Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.**
 - **Never peel off the tube.**
- Doing so may cause overheating, ignition, or explosion.

WARNING

After removing the rechargeable battery, do not keep it within reach of children and infants.

- The battery will harm the body if accidentally swallowed.

If this happens, consult a physician immediately.

If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.

- Failure to do so may cause eye injury.

Consult a physician immediately.

CAUTION

If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.


- Failure to do so may cause inflammation.

► Notes

- Charge the appliance correctly according to these Operating Instructions. (See page 13 “Charging the toothbrush”.)
- This appliance contains battery that is only replaceable by skilled persons. Please contact an authorized service centre for the details of repair.

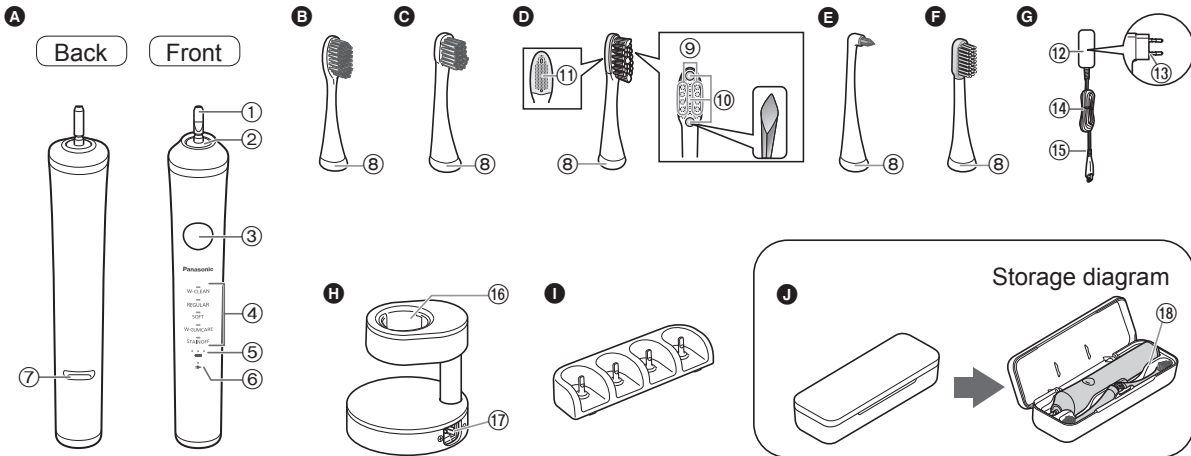Parts identification

A Main unit

- ① Toothbrush shaft
- ② Waterproof rubber seal
- ③ Power switch
(power off/on/mode selection)
- ④ Mode indicator
- ⑤ Battery capacity lamp
()
- ⑥ Charge status lamp
()


- ⑦ Charging ring stopper
(Insert the main unit into the charging stand up to this position while charging.)
- B Extra fine brush**
⑧ Identification ring
- C Compact multi-fit brush**
- D Stain care brush**
⑨ Extra fine side bristles
⑩ High-density bristles
⑪ Tongue cleaner

- E Point brush**
- F Silicone brush**
- G AC adaptor (RE7-87)**
(The shape of the power plug differs depending on the area.)
⑫ Adaptor
⑬ Power plug
⑭ Cord
⑮ Appliance plug

- H Charging stand (RC8-3)**
⑯ Charging ring
⑰ Stand socket
- Accessories
- I Brush stand**
- J Carrying case**
⑱ Inner frame

If the top case has come off

1. Align the mounting section of the top case and the bottom case (2 locations) as shown in the figure.
2. Press the top case in the direction indicated by arrows, taking care not to pinch your fingers.

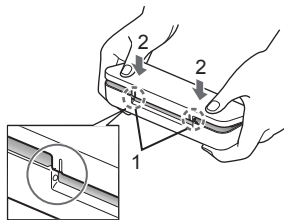Before cleaning your teeth

► This product

The Sonic Vibration Toothbrush is an electric toothbrush with sonic vibration* of approximately 31 000 horizontal and 12 000 vertical brush strokes/minute.

* Sonic vibration in acoustic wave field

► Bleeding of the gums

When you first begin using the toothbrush, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the toothbrush for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as periodontal disease. In this case, stop using the toothbrush and consult your dentist.

Features

► About soft start function

To get accustomed to the high speed vibration, the toothbrush will gradually reach the selected speed and vibration approximately 2 seconds after it starts operating.

► About the 30-second brushing notification timer function

The toothbrush will pause once at 30-second intervals to indicate the brushing time. It pauses twice at 2-minute intervals then starts again.

- This function is not available in W-GUMCARE mode and STAINOFF mode.

► About quick charge

Even if the charge runs out, charging for 2 minutes will provide approximately 2 minutes of use.

- Use time may vary depending on ambient temperature and use conditions.

► About mode storage function

The mode used last is stored.

- The toothbrush operates in the mode last used when the power switch is turned on.


Charging the toothbrush

Preparation

Place the charging stand on a horizontal and stable place.


1 Insert the appliance plug firmly into the stand socket.

- Align with the ⊕ ⊖ symbol as illustrated.

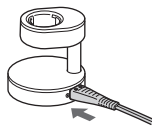
2 With the main unit facing forwards, set straight into the charging ring.

- The charging ring stopper (a) touches the charging ring and the main unit will stop in a raised position.

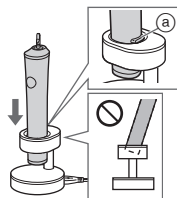► About the approximate remaining charge and charging time

- When you turn off the main unit while the battery capacity is low, the unit beeps.

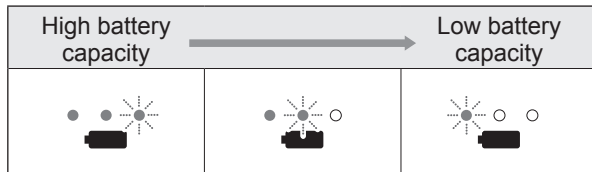
(●:Glow/★:Blinks/○:Turns off)

3 Insert the adaptor into the household outlet.

- Charging is completed after approx. 1 hour.

(The charging time may differ with the charging capacity.)

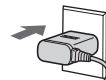► Display while charging

| While charging (Low battery capacity) | While charging | After charging is completed | 5 seconds later after charging is completed |
|---------------------------------------|----------------|-----------------------------|---|
| | | | |

(●:Glow/★:Blinks/○:Turns off)

► How to check if charging is completed

Removing and attaching the main unit after charging will cause the charge status lamp (→) to glow and then turn off after 5 seconds.

4 Disconnect the adaptor from the household outlet.

(for safety and to reduce energy usage)

- Charging the battery for more than 1 hour will not damage the main unit.

► Battery runtime on a full charge

| Mode | Battery runtime |
|-----------------------------|--------------------|
| W-CLEAN W-GUMCARE | Approx. 45 minutes |
| REGULAR SOFT STAINOFF | Approx. 90 minutes |

(Use time may vary depending on ambient temperature and use conditions.)

Notes

- If the charge status lamp (➡) blinks quickly (twice per second), it may be due to a charging error or because the main unit is not properly mounted on the charging stand. Confirm the recommended temperature and mount the main unit on the charging stand properly.
- Recommended ambient temperature for charging is 5 °C – 35 °C. The battery may not charge properly or not at all under extreme low or high temperatures.
- The charging time may vary when using for the first time, or when it has not been in use for more than 6 months. Also, the charge status lamp (➡) may not glow for a few minutes after charging starts, but it will glow if charging continues.


Basic usage

Preparation

Select the appropriate brush for the desired use. (See page 16.)

1 Attach the brush head to the main unit.


- Align the front of the main unit and the front of the brush head.
1. Ensure that the direction is correct.
 - Do not attach the brush head by twisting it into place. Doing so may damage the brush head and main unit.
 - The toothbrush shaft will move when a force is applied. This is required for operation, and is not a malfunction.

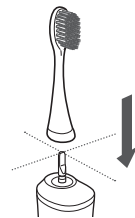2. Apply a small amount of toothpaste.

There is a small gap between the brush head and the main unit (a few mm), but this is not a malfunction.

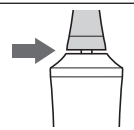
2 To avoid spattering of toothpaste, put the brush head in your mouth and press the power switch to turn on the main unit.

► Select the desired mode

- The device starts in the previous mode used when the power is turned on.
- Press the switch within 5 seconds to change modes. (The blue mode indicator illuminates to indicate the selected mode.)

To switch the mode, press the power switch within 5 seconds.

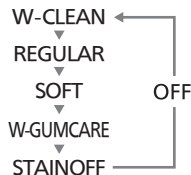
3 Lightly press and move the toothbrush against your teeth.

► Example brushing sequence

Divide your teeth into 4 parts (upper and lower, left and right), and brush one part in order such as “① front of the teeth”, “② biting plane of the teeth”, and “③ back of the teeth” spending 30 seconds. Repeat this same procedure for the other 3 areas.

- By brushing the above 4 parts for 30 seconds each (about 2 minutes), you can evenly brush your teeth without missing any spots in your mouth.
- See page 12 about “30-second brushing notification timer function.”

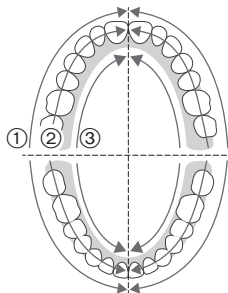4 Turn off the main unit before removing it from your mouth to prevent saliva from spattering.

- Press the switch after at least 5 seconds of operation in one mode to turn the toothbrush off. If you press the switch before 5 seconds of operation have elapsed, the mode changes. (See this page.)

Notes

- Do not put the main unit in contact with hair products, hand creams, etc. Doing so may damage the exterior.
- The appliance may not operate in an ambient temperature approximately 5 °C or lower.
- Use REGULAR mode if the brushing feels too strong in W-CLEAN mode, and use SOFT mode if it feels too strong in REGULAR mode.
- When first using this appliance, some people may experience a tickling or tingling sensation caused by the sonic vibrations. The sensation should subside after you use the toothbrush several more times.
- If you push the brush head too hard, the sound may change, but this is not an abnormality.

Tips on usage

Change the orientation and manner in which you brush your teeth in accordance with your objective and the area you are brushing.

Periodontal pocket care

Recommended mode: REGULAR

Extra fine brush

When you want to brush mainly the periodontal pockets or between the teeth and gums.

- You can also brush the teeth front surface and the biting surface.

Dental plaque removal

Recommended mode: W-CLEAN

Compact multi-fit brush

When you want to make your teeth front surface and biting surface clean.

- You cannot brush the periodontal pocket.

Stain care

Recommended mode: STAINOFF

Stain care brush

When you want stain care during normal brushing.
(This does not mean that all stains can be removed.)

- The intensity of vibration in STAINOFF mode reduces approximately every 30 seconds and then gradually increases.

How to apply/move the toothbrush

Apply the bristles of the brush lightly against the teeth and move the toothbrush slowly.

Front surface and the biting surface of teeth

- Apply the brush head at a 90° angle.

Periodontal pocket and the boundary between the teeth and gums

- Use the extra fine brush.
- Place the brush head at a 45° angle.
- Apply the bristles lightly into the periodontal pocket and move the toothbrush slowly.

For the back of the front teeth or misaligned teeth

- Place the brush vertically.
- Move toothbrush to match the angle of the teeth.

Gum care

Recommended mode: W-GUMCARE

Silicone brush

- When you want to clean the gums.
- Do not use toothpaste.

How to apply/move the toothbrush

Use the toothbrush so that the tips of the silicone brush bristles are lightly touching your gums. Move the brush back and forth gently across your entire gums.

Gum

- Use as close as possible to a right angle.
- Use gently at points of concern. (Use for around 10 to 20 seconds at each area.)
- Recommended usage is once a day, for approximately 2 minutes.

Notes

- Be careful when using the brush head. Forcefully pulling or biting it may tear it.
- Do not apply the brush to your teeth. (This can cause the silicone brush to degrade.)

Tongue care

Recommended mode: SOFT

Tongue cleaner

Stain care brush (Back)

When you want to clean the tongue.

How to apply/move the toothbrush

Apply the tongue cleaner so that it barely touches the tongue, and move forward with slight force from the area that becomes dirty.

Tongue

- Repeat several times.
- Poke your tongue out forward and apply the brush while you are breathing out to prevent emetic response.

Note

- You may injure your tongue if the tongue cleaner is pressed against the tongue too strong.

Point brushing

Recommended mode: W-CLEAN

Point brush


For the gap between teeth, the rear of the back teeth or misaligned teeth.

How to apply/move the toothbrush

Apply the bristles of the brush lightly against your teeth.


Gap between teeth

- Place the brush at gap between teeth or areas that are hard to brush.
- Move the toothbrush slowly in line with the gums.

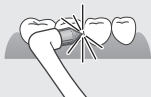
Rear of the back teeth

- Apply the brush to the rear of the back teeth thoroughly, as they are hard to reach.

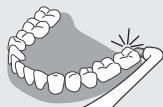
Misaligned teeth

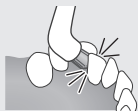Braces

- Apply the brush into the gap between braces and teeth.

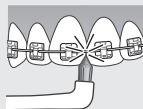
After use

After using the toothbrush, clean the toothbrush to maintain cleanliness of the device.

- Use neutral kitchen detergent when cleaning the toothbrush with detergent.
- Do not use thinner, benzine, alcohol, etc. to clean the appliance.
Doing so may cause failure, cracking or discoloration of the parts.

Remove the brush head.

Remove the brush head straight from the toothbrush shaft.

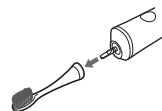Cleaning the brush head

Wash the brush head under running water.

Wash the soiling on the inside and bottom areas of the brush head thoroughly.


- Do not forcefully pull the bristles of the brush or bend them out excessively.
- Do not wash the tongue cleaner excessively with pointed objects or your fingernails.
- Do not wash with water hotter than 80 °C.
Doing so may damage the brush head.

Cleaning the main unit

Wash the main unit under running water.


- Do not damage or rub the waterproof rubber seal (A) with hard materials.
Do not leave the main unit under water. Water may enter inside.
- Do not wash with water hotter than 40 °C.
Doing so may damage the main unit.

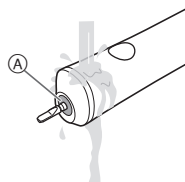Cleaning the AC adaptor and the charging stand

Remove dirt with a dry cloth.

- Do not wash with water.

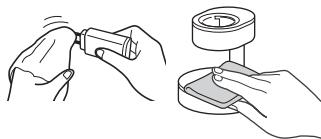Cleaning the brush stand


Wash and wipe off water.


Cleaning the carrying case

1. Remove the inner frame.

- ① Hold the right end of the carrying case.
- ② With the right end of the inner frame facing up, pull up the inner frame.

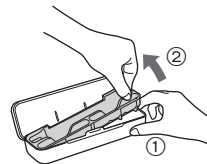2. Wash with water and wipe off water.

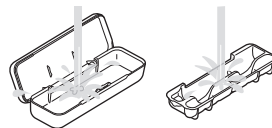
► After cleaning

Wipe off water from the brush head and main unit for storage.

- Do not charge while it is wet. Doing so may cause soiling.

► If the identification ring comes off

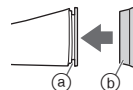Attach by aligning the indentation (a) in the brush head with the raised part (b) of the inside of the identification ring.

- Do not forcefully pull the identification ring.
Doing so may cause damage.
- Using the toothbrush without the identification ring has no effect on its performance.

Changing the brush head (Sold separately)

- Brushes are consumables.
- The brush should be replaced with a new one every 3 months or so due to hygiene reasons, even if its shape is unchanged.
Even within 3 months, please exchange the brush head if the bristles are frayed.
- Colored identification rings are included with the separately sold brush.
Optional brush heads with identification ring to allow several persons to share the main unit.

| Part number | Replacement Brush for EW-DP52 | |
|-------------|-------------------------------|---|
| EW0917 | Extra fine brush |  |
| WEW0915 | Compact multi-fit brush |  |
| WEW0929 | Stain care brush |  |
| EW0940 | Point brush |  |

► Silicone brush

Silicone brush is treated as service parts.
If you want to replace the silicone brush,
please contact an authorized service centre.

FAQ

| | Question | Answer |
|-------------|---|--|
| While using | Why does the toothbrush pause once, or the intensity of vibration reduce, at 30-second intervals? | ► The toothbrush will notify you of the elapsed brushing time every 30 seconds. (See page 12 about "30-second brushing notification timer function.") (See page 16 about "STAINOFF mode.") |
| Battery | Will the battery deteriorate if not used for an extended period of time? | ► The battery will deteriorate if not used for over 6 months, so make sure to perform a full charge at least once every 6 months. |
| | Can I charge the appliance after every use? | ► The battery used is a lithium-ion, so charging after every use will have no effect on battery life. |

| | Problem | Possible cause | Action |
|----------------|--|---|---|
| While charging | The charge status lamp (➡) does not glow when charging starts. | <ul style="list-style-type: none"> ▶ You have just purchased the toothbrush or have not used it for a while. ▶ The main unit is not properly mounted on the charging stand. | <ul style="list-style-type: none"> ▶ The charge status lamp (➡) may not glow for several minutes after starting to charge it, but it will glow if you continue to charge it. If the charge status lamp (➡) does not glow after a while, consult with an authorized service centre for repair. ▶ Insert the main unit in the charging ring and check that the charge status lamp (➡) glows. (See page 13.) |
| | The charge status lamp (➡) is blinking rapidly. (Twice every second) | <ul style="list-style-type: none"> ▶ The temperature is too low or high while charging. | <ul style="list-style-type: none"> ▶ Charge the toothbrush within the recommended charging temperature of 5 °C to 35 °C. ▶ If the charge status lamp (➡) blinks even when charging within the recommended temperature range, please contact an authorized service centre for repair. |

Battery life

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge.

(The operating life of the battery may differ depending on usage or storage conditions.)

Removing the built-in rechargeable battery

Remove the built-in rechargeable battery before disposing the toothbrush.


The battery is to be disposed of safely.

Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing the toothbrush, and must not be used to repair it. If you dismantle the toothbrush yourself, it will no longer be waterproof, which may cause it to malfunction.

- Press the power switch to turn on the power and then keep the power on until the battery is completely discharged.

- Follow the procedure below, using a screwdriver to disassemble the appliance. Perform steps ① to ④ and lift the battery, and then remove it.
- Take care not to short circuit the positive and negative terminals of the removed battery, and insulate the terminals by applying tape to them.

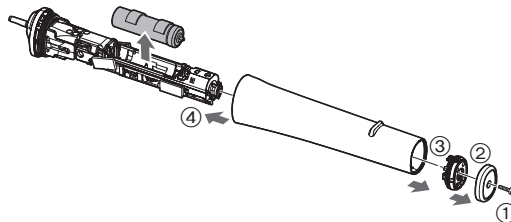For environmental protection and recycling of materials

This toothbrush contains a Li-ion battery.

Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

Specifications

| | |
|---------------------------|---|
| Power supply | See the name plate on the product (Automatic voltage conversion) |
| Charging time | Approx. 1 hour (Room temperature: 5 °C – 35 °C) |
| Battery | 1 Lithium-ion battery (3.6 V; 520 mAh) |
| Airborne Acoustical Noise | 59 (dB (A) re 1 pW) |

This product is intended for household use only.

For more information about the energy efficiency of the product, please visit our website, www.panasonic.com, and enter the model number in the search box.

Disposal of Old Equipment and Batteries Only for European Union and countries with recycling systems

These symbols on the products, packaging, and/or accompanying documents mean that used electrical and electronic products and batteries must not be mixed with general household waste.

For proper treatment, recovery and recycling of old products and batteries, please take them to applicable collection points in accordance with your national legislation. By disposing of them correctly, you will help to save valuable resources and prevent any potential negative effects on human health and the environment.

For more information about collection and recycling, please contact your local authority. Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

Note for the battery symbol (bottom symbol)

This symbol might be used in combination with a chemical symbol. In this case it complies with the requirement set by the Directive for the chemical involved.

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© Panasonic Corporation 2019-2022
Panasonic Marketing Europe GmbH
Panasonic Corporation
AB Yetkili Temsilcisi: