
# A E N O

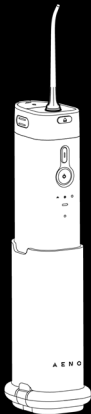Quick Start Guide

# Dental Irrigator D12


(Model: ADI0002)


v 2.1.0 | 02.04.2024


ENG.....	14
ARA.....	18
BOS.....	21
BUL.....	25
CES.....	29
DEU.....	33
ELL.....	37
EST.....	42
FRA.....	46
HRV.....	50
HUN.....	54
HYE.....	58
ITA.....	62
KAT.....	66
KAZ.....	70
LAV.....	74
LIT.....	78
NLD.....	82
POL.....	86
POR.....	90
RON.....	94
RUS.....	98
SLK.....	102
SLV.....	106
SPA.....	110
SRP.....	114
UKR.....	118
UZB.....	122

# A

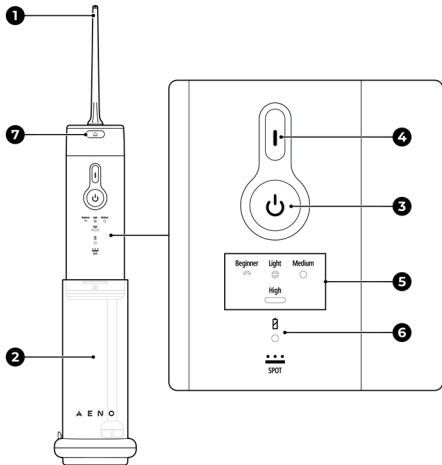
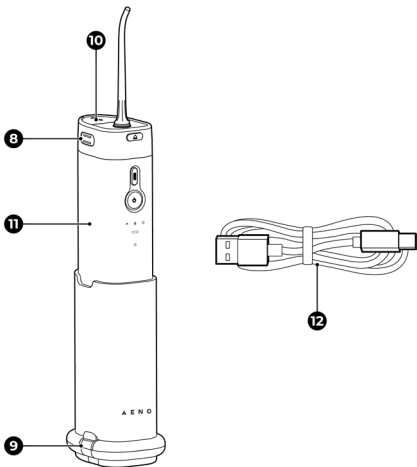
**C**

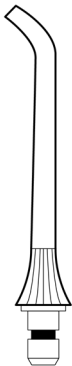
**13**

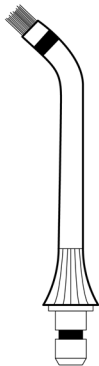
**14**

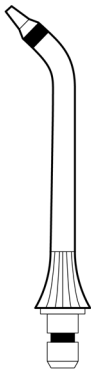**15**

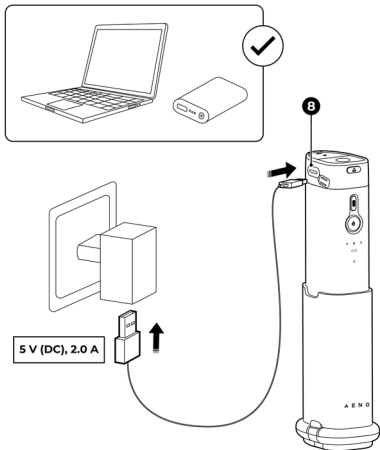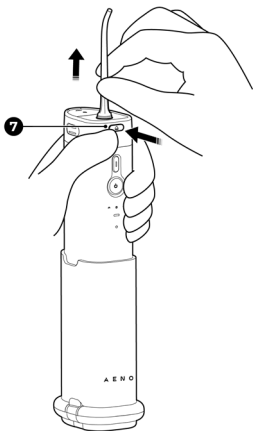
**E****F**



**G**



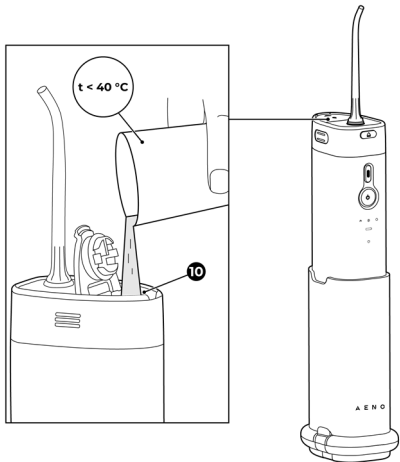
**H**

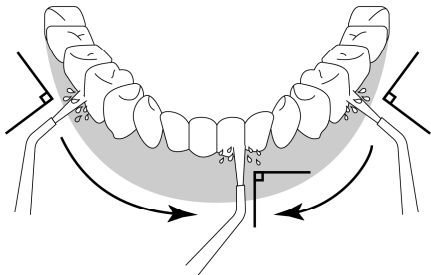I

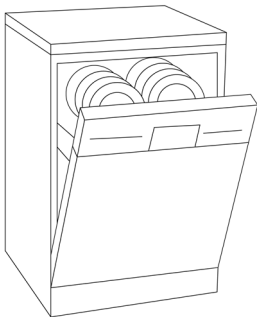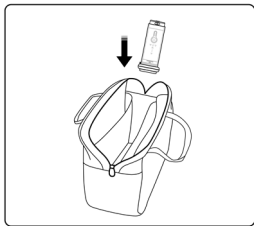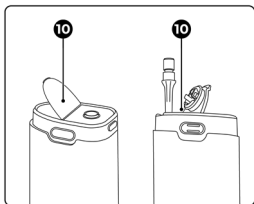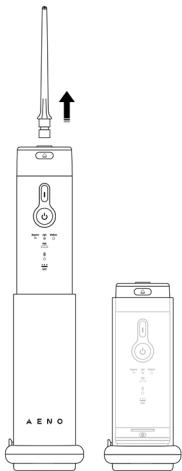**ENG** **AENO™ ADI0002** Dental Irrigator is designed for oral care and removing plaque and food debris from hard-to-reach areas.

### **Technical Specifications**

Power supply (input): 5.0 V (DC), 2.0 A. Power: 10.0 W (max.). Battery: non-removable, lithium-ion, 2,500 mAh, 3.7 V. Noise level (at a distance of 0.3 m): ≤ 72 dB. Water pressure: 206–965 kPa. Pulsation: 1,900 impulses/min. Full battery charge time: 3 h. Water tank: 200 ml. Ingress protection rating: IPX7. Operating time on a full charge: at least 200 min. Dimensions when folded (L×W×H): 71×52×131 mm. Operating conditions: temperature +5...+40 °C, relative humidity ≤ 75 %. Storage conditions: temperature -10...+60 °C, relative humidity ≤ 75 %.

### **Scope of Supply** (see Fig. **B**, **C**)

**AENO™ ADI0002** Dental Irrigator (White) (11), 3 replaceable nozzles (13–15), USB (A) / USB (C) cable (0.6 m) (12), Quick Start Guide.

### **Limitations and Warnings**

The device is not intended for use by persons under the age of 6. Supervise the use of the device by children, as well as by persons with physical or mental disabilities. Do not allow children to play with the device. If you have an oral disease or a pacemaker, consult a doctor before using the device. Do not use the nozzles or device if they are damaged. Do not leave the device that has been plugged in unattended. Use the irrigator, as well as clean it and its accessories only after disconnecting the device from the power supply. **Use purified water only without any mouthwashes and other liquids or substances.** Do not point the outlet of the nozzle into your eyes while the device is operating. Replace nozzles at least once every 3 months.

If you have any issues or difficulties using your AENO™ device, please email support at [support@aeno.com](mailto:support@aeno.com) or chat with us at [aeno.com/service-and-warranty](https://aeno.com/service-and-warranty). Our specialists will help you solve the problem, and you will not have to spend time and effort visiting the store.

## Preparing the Device for Operation

Unpack the device. Wipe its casing with a cloth dampened with water, then wipe dry. Detach the water tank (see Fig. **F**). Rinse and dry the water tank and all the nozzles as described in the “Cleaning and Maintenance” section (“Daily care”).

## Charging the Device

**Note.** The power adapter is not included in the package. You can use any power supply (adapter, portable charger, etc.) as long as it ensures 5.0 V (DC), 2.0 A power parameters.

Connect the device to a power supply as shown in Figure **D**. Statures of the battery charge indicator (6) while charging:

- The device is charging: the indicator is blinking;
- The device is fully charged: the indicator is lit up steadily.

When charging is complete, disconnect the device from the power supply, detach the cable (12) and close the USB port (8) by firmly inserting the cover.

Statures of the battery charge indicator (6) while in operation:

- Less than 10 % – the indicator is blinking, the device does not turn on;
- 10–100 % – the indicator is lit up steadily.

Attach one of the replaceable nozzles (13–15) to the casing by pushing it until it clicks into place. The nozzle can rotate 360°.

**Intended Purpose of the Irrigator Nozzles** (see Fig. **C**):

(13) – standard nozzle for removing food debris and plaque;

(14) – orthodontic nozzle for braces and dentures;

(15) – nozzle for cleaning periodontal pockets.

Fill the water tank (2) as shown in Figure **G**. Fasten the cap of the water filling opening (10). Afterwards, select the operating mode of the device.

## Selecting the Operating Mode of the Device

**Note.** To select the operating mode, point the outlet of the nozzle in a safe direction so that the jet of water does not damage your surroundings and you can see the mode selection panel (5).

Press the button (3) to turn on the device. The operating mode indicator on the irrigator panel (5) will light up. Use the mode selection button (4) to select the desired mode.

### **Operating Modes of the Device**

- “Beginner” – gentle cleaning for sensitive teeth and gums. Also recommended for those using the device for the first time;
- “Light” – uniform cleaning before moving on to intensive cleaning;
- “Medium” – daily standard cleaning, plaque removal;
- “High” – intensive cleaning, thorough oral care;
- “Spot” – removal of food debris from interdental spaces, prevention of tartar formation.

**Note.** The “Spot” mode does not have its own indicator on the mode selection panel (5). When the device is operating in this mode, all four panel indicators are lit up.

By default, the device is set to “High” operating mode the first time it is turned on. Subsequent single presses of the mode selection button (4) alternately change the operating modes, which are displayed on the panel (5). The last selected mode is saved in the device memory and will be automatically selected the next time the irrigator is turned on. After selecting the mode, move on to device operation.

### **Device Operation**

Direct the outlet of the nozzle (1) into your mouth cavity. Tap button (3) to turn on the device. Move the outlet of the nozzle (1) along the gum contour as shown in Figure H. During cleaning, you can switch the operating modes using button (4). After cleaning is complete, tap button (3) to turn the device off. If the device is not turned off after use, it will automatically switch off after no more than 2 minutes. After each use, rinse the device and the used nozzles as described in the “Cleaning and Maintenance” section (“Daily care”).

**Note.** It is not recommended to use nozzle (15) in the “Spot” mode to avoid injuring periodontal pockets.

## **Cleaning and Maintenance**

- Daily care. Open the cap of the water outlet (9) and drain the remaining water from the tank. Detach the nozzle using button (7) as shown in Figure **E**. Rinse the tank and nozzles under warm running water without using any cleaning agents and dry them thoroughly.
- Total care. Perform the total care of the water tank once per week. To do this, detach the tank as shown in Figure **F**. Rinse its inside and outside under warm running water using a soft brush or sponge, without cleaning agents, and dry it thoroughly.

## **Troubleshooting**

- The device does not turn on. Possible reason: the battery is low. Solution: charge the device for at least 3 hours.
- The device does not eject water. Possible causes: the outlet of the nozzle is clogged; there is no water left in the tank. Solution: clean the outlet of the nozzle; fill the tank with water.

**ATTENTION!** If none of the suggested solutions has helped you solve the problem, contact the supplier or service centre. Do not disassemble the device or attempt to repair it yourself.

ASBISC reserves the right to modify the device and to make changes and additions to this document without prior notice to users.

The warranty period and service life is 2 years from the date of retail sale of the device. Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Made in China. All trademarks and brands mentioned herein are the property of their respective owners. The date of manufacture is indicated on the package.

Up-to-date information and a detailed description of the device, as well as connection instructions, certificates, information on companies that accept quality claims and warranties are available at **[aeno.com/documents](http://aeno.com/documents)**.

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**ENG Recycling information** These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or turned in to a local recycling facility. You should contact your local household waste disposal service for details.

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