MUSCLE MASSAGE GUN

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The Muscle Massage Gun PG-M190 is designed with 3 working modes for more relaxed service in moderate strength; in addition, it is provided with several massage heads for your free choice and different service effects; its built-in large-capacity rechargeable batteries provide it a long-term service life and your use anytime. For its delicate and fashionable design, and convenience in carrying, it will be a good assistant for your relaxation at home or in travel.


PRODUCT FEATURES

- Simulation of manual beating to deeply massage your fascia and relax your muscles.
- Mini-type design and low weight for easy carrying.
- Able to solve the problem of getting tired due to long-time holding of the thick handle of traditional fascia massager.
- Massage muscle tissues to relieve fatigue.
- Reasonable 6mm massage stroke to provide a comfortable experience.
- 3 automatic working modes for free choice according to your needs and being suitable for a wider range of users.
- 4 massage heads to ensure each part massaged, so as to meet different needs in daily relaxation; brushless tailor-made motor to provide surging power with low noise.
- Built-in rechargeable batteries for power supply and recharging with USB5.0V.
- Simple, elegant and fashionable appearance to make it easy to carry and massage anytime anywhere.

PRODUCT FEATURES

- 1 main unit (including 1 ball-type massage head)
- 3 non-spherical massage heads
- 1 USB cable
- 1 copy of user manual

PRODUCT DIAGRAM

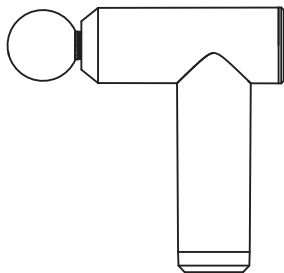Main unit

Ball-type massage head

Suitable for massaging large muscle groups such as arms, back, thighs, and calves.

U-type massage head

Suitable for massaging both sides of the neck and spine, and also able to be used on shoulders, triceps, etc.

Bullet-type massage head

Suitable for massaging deep tissues such as meridians, joints, palms, and soles.

Flat massage head

Suitable for relaxing and remodeling muscles.

HOW TO OPERATE IT


How to install the massage head:

Select a massage head based on your needs to massage different body parts, and then insert it into the round hole on this product, as shown in Figure 1; Insert it completely into the round hole to complete the installation, as shown in Figure 2.


1. When it is off, long press "On/Off" to turn on it; the massage strength defaults to Level 1.
2. When it is on, long press "On/Off" to turn off it.
3. When it is on, short press "On/Off" to change the massage strength (Level 1 – Level 2 – Level 3 – Level...in cycle).
4. When it is switched to Level 1, the Level 1 indicator will be always on; when it is switched to Level 2, the Level 2 indicator will be always on; when it is switched to Level 3, the Level 3 indicator will be always on.
5. It will stop operating and fail to be used in charging.

COMMONLY MASSAGED MUSCLES AND RECOMMENDED MASSAGE DURATIONS AND TIMES

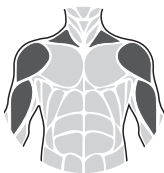
Shoulder muscles
15~20s x 3 times

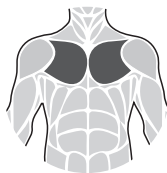Pectoral muscles
30s x 3 times

Back muscles
30s x 3 times

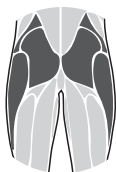Arm muscles
30s x 4 times

Thigh muscles
30s x 3 times

Calf muscles
30s x 3 times

TIPS

1. Reasonable use of this product can relieve muscle fatigue, but the actual effect will vary from person to person depending on their body conditions;
2. Never use this product on the same body part for more than 30 seconds;
3. This product has no effect of treating diseases for it not being a medical device.

CLEANING AND MAINTENANCE

1. Be sure to turn off this product before cleaning, and please clean it when it is not in charging.
2. Please use neutral detergent and a soft cloth for daily cleaning, with any corrosive detergent prohibited; please never soak or rinse it in water.
3. Water and other liquids should be prevented from penetration into the main unit of this product during cleaning, so as to avoid the risk of short circuit and electric shock.
4. Please place this product in a safe place after use, and do not place any heavy object on it to prevent damage or deformation.
5. Please keep this product out of the reach of children.
6. Please store it in a dry and well ventilated place to avoid strong sunlight and corrosion by corrosive gases or liquids.
7. Please power off this product when not in use for a long time, and charge it once every two months.
8. Please keep this product and its power cord separated when not in use.

TROUBLESHOOTING

Please don't disassemble or attempt to repair this product by yourself.


Trouble	Inspection item	Solution
Failure to star	The built-in rechargeable battery is in a low level?	If so, please charge it with the USB cable.
Failure to be charged or no charging display	The USB cable is properly connected?	Check both ends of the USB cable to see if it is properly connected.
	The charging source is in a low voltage?	Check if the charging source is in normal power output.
Abnormal vibration	Is the massage head properly installed?	Please install the massage head again.

USE OF THE RECHARGEABLE BATTERY

1. This product is designed with a built-in rechargeable battery. Please charge it timely when it is in a low battery.
2. Please use an adapter with a power output of 5V---1A to charge this product. To charge it, please insert the USB end of the charging cable supplied into the adapter, and the other end into its Type-C charging port.
3. Charging indication: the red light flashes during charging and stays on after fully charged.

WASTE DISPOSAL

Please subject to relevant urban regulations on environmental protection to dispose this product, otherwise the environmental pollution may be caused.

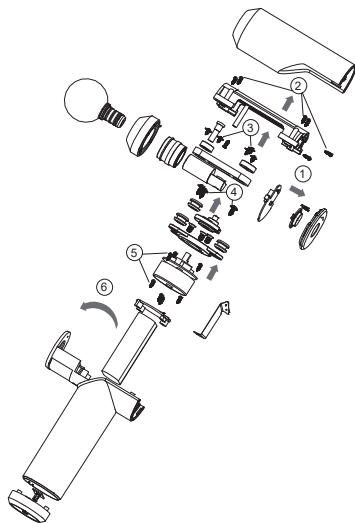ATTENTIONS

1. Please read and understand all the contents of the User Manual before use of this product.
2. Those who do not have the consciousness or ability to express, such as children, are prohibited from use of this product.
3. In case of any discomfort or skin abnormality in use of this product, please stop using it immediately and then consult with your doctor.
4. Before use of this product, please check first if the massage head is properly installed. The massage effect will be affected with a massage head not properly installed.
5. During normal use of this product, please never press the massage head forcefully.
6. To avoid damages, please never aim the massage head at human joints during use of this product.
7. Please never have this product used by those with physical, sensory or mental disabilities, or lack of relevant experience and knowledge (including children).
8. Please never use this product to body parts with bruise, pain or swelling.
9. Please never use this product to body parts such as your head, neck, thoracic cavity, and armpit where the muscles are thin and close to organs, vital tissues, or the aorta.
10. During use of this product, please be careful not to get it close to your hair to avoid your hair being caught by the massage head.
11. Please don't use this product to the same body part for a long time. It is recommended to control the massage time within 30 seconds. After 30 seconds, you should pause for a while before continuing the massage, or use it to other body parts.
12. For each use of this product, please never make it operating for more than 10 minutes, so as to avoid it to heat up due to excessive use.
13. This product is prohibited from being immersed or placed into water for cleaning.

14. Please never put this product into fire or place it near fire.
15. It is prohibited from using this product in high temperature, flammable and humid places or those with electromagnetic radiation.
16. It is forbidden to disassemble this product by yourself to avoid failures caused due to collision or water penetration.
17. When not in use for a long time, please try to charge it once a month to avoid battery aging.
18. Please never place this product on an unstable table to avoid it being dropped and damaged.
19. Please keep the plastic bag used to package this product out of the reach of children to avoid suffocation and other dangers caused by playing with it.

WARNINGS

1. Pregnant women, infants and young children are prohibited from using this product; and menstruating women are not recommended to use it.
2. Patients with tumors or cancer are prohibited from use of this product.
3. Patients after a surgery (including heart bypass, stent implantation, etc.) are prohibited from use of this product.
4. Patients sensitive to hyperthermia, or with skin sensory disturbance or abnormalities are forbidden to use this product.
5. Patients with a moderate and severe cardiovascular or heart disease are forbidden to use this product.
6. For patients with edema caused by acute trauma, bacterial infection, or other diseases, please consult a doctor before use of this product.
7. This product cannot be used for diagnosis and treatment of diseases for it not being a medical device.

PRODUCT FEATURES

- Simulation of manual beating to deeply massage your fascia and relax your muscles.
- Mini-type design and low weight for easy carrying.
- Able to solve the problem of getting tired due to long-time holding of the thick handle of traditional fascia massager.
- Massage muscle tissues to relieve fatigue.
- Reasonable 6mm massage stroka to provide a comfortable experience.
- 3 automatic working modes for free choice according to your needs and being suitable for a wider range of users.
- 4 massage heads to ensure each part massaged, so as to meet different needs in daily relaxation.
- Brushless tai lor-made motor to provide surging power with low noise.
- Built-in rechargeable batteries for power supply and recharging with USB 5. 0V.
- Simple, elegant and fashionable appearance to make it easy to carry and massage anytime anywhere.

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