Sonic Toothbrush Family

Sonic Toothbrush Family White/Black

Meriden Sp. z o.o.
ul. Piastów96
42-290 Blachownia, Poland
NIP: 9492217625
meriden@meriden.pl
www.meridenglobal.com

www.facebook.com/meridenglobal
www.instagram.com/meridenglobal

Designed by Meriden in Poland
Made in P.R.C.

Model: MD 400W, MD 401B

Before using this product,
please read carefully
this user's manual
and keep it properly
for later reference.

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MD 400W, MD 401B

Information about a Sonic Toothbrush FAMILY WHITE or BLACK

All-in-one cleaning and protection system

Product type	MD-400W/ MD-401B
Working mode	Sensitive/Clean/White
Working voltage	4V ---
Rated power	3W
Button type	Touch button
Display mode	LED light
Operating time	240min
Charging time	15hrs
Charging current	55mA
Battery specification	Lithium battery/14500/600mAh
Power reminder	Indicator light flashing at low power
Waterproof level	IPX7
Executive Standards	GB4706.1-2005 GB4706.59-2008 Q/RMZN 1-2016

Please, be informed that brush head parts are made of food grade material.

Brush head

Brush head cover

Toothbrush holder

Touch button

Long press to turn
on / off

Mode switch

Clean mode

White mode

Sensitive Mode

Power

Induction
charging base

Start


EN

Five steps to protect your mouth

EN


01 Installing

Align toothbrush body with brush head and install

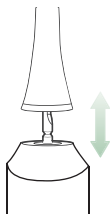02 Rinsing

Rinse with water and squeeze appropriate amount of toothpaste.

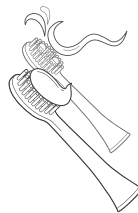03 Selection

Personalized teeth protection with 3 cleaning modes

Tip: Select Sensitive mode when you first use the toothbrush.

Sensitive

62 000 reps/min
Gentle care

Clean

82 000 reps/min
Cleaning and protection


White

96 000 reps/min
Tooth wash level
- White

04 Brushing


Keep bristles against your gum at a 45° angle and clean all four areas

Tip: Put the brush head into your mouth before turning on.

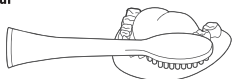05 Brushing

Gargle with water and rinse the brush head.

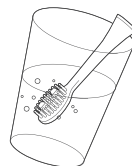Smart

EN

Know your teeth better

5 seconds of smart memory

Smart and considerate in remembering your habits.

2 minutes automatic timer

Scientifically plan your brushing time.

Reminds to change area every 30 seconds

Clean every corner of the mouth in a balanced manner.

Bristles wear out after 3 months.

Change brush heads regularly.

60 days of exceptional endurance

With one charging the toothbrush may be used for a period of two months.

1-start mode change

Easy and simple One-key control.

Tips

EN

Habits that protect your mouth

01 Deep cleaning

Brush your teeth every morning and evening for at least 2 minutes each time.

02 Brushing at the right time

Brush your teeth 20 minutes after a meal to protect your enamel.

03 Gentle gums

Use soft bristles to clean sensitive gums.

Dental health recommendations

Children

Start as a good example for children. Help them to protect their teeth as well as gums and reduce incidence of dental caries.

- Brush your teeth every morning not less than 2 minutes each time.
- Ensure a balanced diet and control snacking.
- Use toothpaste to prevent dental caries.
- Use toothbrush with special bristles for children to gently clean every corner of teeth.
- After brushing your teeth use the toothbrush bristles to brush your tongue to effectively reduce plaque.
- Regular oral examination.

Teenager

Be an owner of neat and beautiful teeth. More and more teenagers choose to have orthodontic appliances which makes cleaning tooth harder. Therefore, daily care cannot be treated lightly.

- Brush your teeth every morning not less than 2 minutes each time.
- Use professional toothpaste to remove plaque and prevent dental caries.
- Gargle your throat after meals.
- Regular dental examination.

The elderly

By maintaining correct dental care habits healthy teeth may last your entire life. It is important to care for your mouth no matter how old you are.

- Constantly use toothpaste to prevent dental caries.
- Drink plenty of water or chew sugar-free gum to increase saliva secretion what reduces risk of dental caries.
- Dentures also require careful and thorough cleaning.
- Gently care for your gums to promote good blood circulation what helps to maintain healthy gums.

Care

EN

Daily care and maintenance

Brush head

Rinse with water and leave to dry after each use. Please, use original replacement brush heads for your own safety and to extend the life of the product.

Toothbrush body

The new Tooth Brush will have residual amount of electricity after leaving the manufacturer. Please use up the rest of electricity before charging. Charging indicator will flash continuously to indicate low power. After using up the electricity charge again to extend the service life of the lithium-ion battery. During charging process the indicator light will flash in a form of breathing light and stay on when charging process is complete. While charging do not place any metal objects on the charging base to avoid danger.

Do not drop

Do not use the toothbrush if toothbrush body or brush head is damaged.

Keep away from heat sources

The toothbrush or charging base cannot be used in high temperature environment.

Do not disassemble

Please contact after-sales department for any repairs.

Q&A

EN

Q&A Notes

Q: Why does toothbrush vibration seem to be getting weaker?

A: Connection between the brush head and the body is too tight, please remove the brush head and then re-install.

Q: I have charged the toothbrush but why is operating time still very short?

A: Please make sure that the toothbrush is placed vertically when charging. Incorrect positioning leads to poor charging effect. If the problem has not been resolved, the battery may have reached the end of its useful life. Please contact the After-sales Department of Meriden for replacement.

Q: When charging an electric toothbrush which was not used for a period of time, indicator light does not shine.

A: When a toothbrush is not used for a period of six months or longer the battery enters a dormant state and needs 6-20 minutes after re-charging to wake up. It is a normal phenomenon.

Q: What to do if electric toothbrush does not work?

A: The toothbrush battery might be exhausted and continue to work normally after charging. If the toothbrush does not work after charging, please contact After-sales Department of Meriden for repair or replacement.

Q: Can the power cord be replaced?

A: The power cord cannot be replaced. If the cord is damaged, the electric toothbrush should be discarded.

Service After-sales service card

Notes: Please fill out the warranty card and keep it for future use.

This product can enjoy two years of free warranty service counting from the date of purchase. Free warranty service does not apply in the following cases:

1. Damage caused by improper use or improper storage;
2. Self-disassembly of the product or when product was repaired outside of company's special point of repair;
3. Lack of valid proof of purchase or the valid proof of purchase was altered;
4. Product damage due to voltage instability or exceeding of the specified range;
5. Failure to operate according to the instruction manual;
6. Product exceeded warranty period or scope of warranty;

Model of the product
Date of purchase
Customer name
Contact information
Fault description

Tips: Warranty covers the toothbrush body and charger (excluding the cover of brush head).

Meriden Sp. z o.o. reserves the right of final explanation.

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