
meriden

Water Flosser Home & Travel

Before using this product, please read carefully this user's manual and keep it properly for later reference.

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
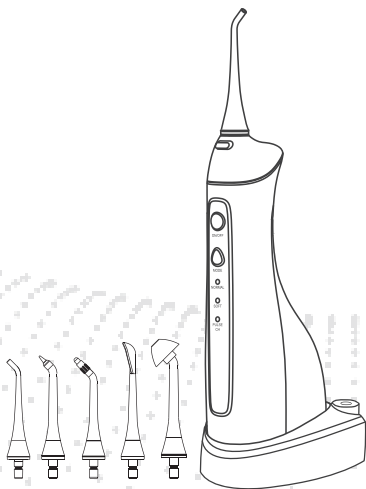
MD 1622

Safety Warning


To prevent burning, electrocution, fire or injury to people, take the following safety precautions:

- Do not plug this device into a voltage system different from the voltage system specified on the device or charger.
- Do not touch the charger with wet hands.
- Do not use the product if there is any damage to the cord, plug or the product itself, or the product is not working properly.
- Do not direct water under the tongue, into an ear, nose or other delicate areas, as the pressure may injure them seriously.
- Use this product only as indicated in these instructions or as recommended by your dentist.
- Fill the reservoir with water or dentist-recommended solutions only.
- Be careful not to drop the product. Never push any objects through the openings in this product. Power cord cannot be replaced- in the case of the power cord damage, the product must not be used.
- Keep the charging unit and cord away from heated surfaces.
- Do not use this product outdoors or operate it in the area where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use iodine, bleach or water insoluble concentrated essential oils in this product. The use of such substances may lead to impairment of the product performance and shortening of the product life.
- Remove any oral jewellery before using this product.
- Do not use this product if you have any open wounds on the tongue or in the mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Children should be supervised to ensure that they do not play with the device.

Product description

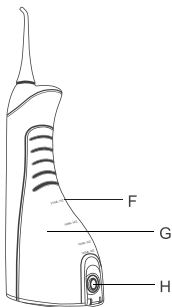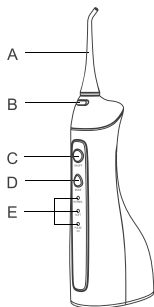- A. Tip
- B. Tip eject button
- C. ON/OFF button
- D. Mode switch
- E. LED indicators
- F. Measuring lines
- G. Reservoir
- H. Reservoir fill door
- I. Inductive charging base

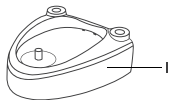
1. Standard jet tip: Polishes teeth and removes waste from between the teeth.
2. Special jet tip for periodontal pockets: Prevents dental calculus and dental plaque.
3. Special jet tip for orthodontics: Cleans hard to reach areas around dental brace.
4. Jet tip for tongue cleaning: Removes the oral malodor and helps to keep fresh breath.
5. Special jet tip for nose cleaning: Reduces rhinitis symptoms and helps to keep nasal cavity healthy.

How to use


Charging the unit

1. Charging indicator light will turn red and blink when the battery is in the need of charging.
2. Plug the power cord into the wall socket (100-240V, 50/60Hz) and place the water flosser on the inductive charging base to charge the unit.
3. While charging, indicator on the handle will flash and when fully charged (approximately 24 hours) it will stay on continuously.
4. A full charge will last more than 30 minutes of standard use.

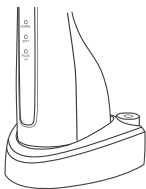Installing and removing the tip


- Insert the tip into the center of the knob at the top of the flosser handle.
- Press firmly until the tip clicks and locks into place.
- To remove the tip from the handle, press the tip eject button and pull the tip from the handle.

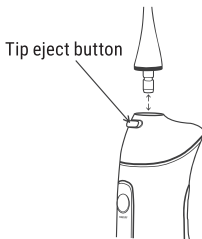Caution: Do not press eject button while the unit is running.


Filling the reservoir

1. Flat on the water flosser keeping the reservoir fill door upwards.
2. Open the reservoir fill door and fill the reservoir with lukewarm water.
3. Close the reservoir fill door when the reservoir has been filled.

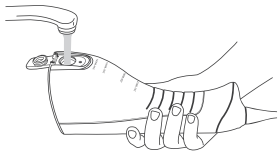Adjusting the pressure settings

- You can adjust the pressure setting on your water flosser by pressing the mode button on the handle.
- Start with the lowest setting (soft<normal<pulse) and increase to your liking.

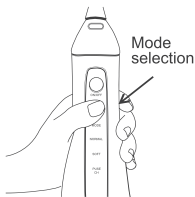
Recommended cleaning tips

- Lean over and place the tip in your mouth. Aim the tip at your gumline and turn the unit ON.
- For best results, start with the back teeth and work towards the front teeth.
- Glide the tip along the gumline and pause briefly between the teeth until you have cleaned the inside and outside of both the upper and lower teeth.
- Direct the jet stream at a 90° angle against your gumline. Slightly close the lips to avoid water splashing out but at the same time allow it to flow freely from your mouth into the bathroom sink.
- For best results, keep the unit in an upright position during use.

After use

1. Turn the unit off.
2. Empty any remaining liquid in the reservoir.
3. Clean the unit directly under the tap as it is IPX7 water-proof design.
4. Dry the unit with a cloth before putting it back on the charging base.

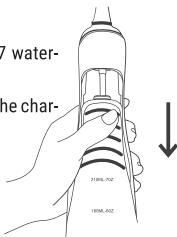Cleaning

- When necessary clean the product with use of a soft cloth and mild non-abrasive cleanser to wipe the product exterior.
- Water reservoir is removable for easy cleaning.
- It is recommended that that water reservoir is removed and cleaned once a week.

Storage

- Store the device out of children's reach.
- Storage space should be dry and clean with no direct sun exposure and away from heat sources.
- If the water flosser needs to be stored for an extended period of time, (over 3 months) be sure to fully charge it.

Troubleshooting

Problem	Possible cause	Solution
Does not work	The product is new - it has just been purchased or it has been idle for more than 3 months	Charge the device for 24 hours.
It only works for a few minutes even after charging	Charging time was insufficient	Charge the unit for at least 24 hours
	The battery has reached the end of its life cycle.	Contact seller for support.
Water pressure is very low	Dirt has stuck in the tip.	Replace the tip
	The tip is deformed	
	Dirt has stuck in the filter	Clean the filter
Water does not come out	The reservoir is empty	Fill up reservoir
	The water irrigator has been tilted too much when used	Use the device in the vertical position

Environmental Protection

If you want to discard of this product, please, do not place it together with domestic waste, but send it to the officially designated recycling centre and dispose of the battery at the officially designated recycling centre. Doing so can support environmental protection.

Warranty and Service

Please contact us by calling 4008002138 if you are need of service, consultation or in case of any problems.

The product comes with 24-month warranty.

Explanation:

Please fill out the warranty card and keep it for future use. This product can enjoy two years of free warranty service counting from the date of purchase. Free warranty service does not apply in the following case of:

1. Damage caused by improper use or improper storage;
2. Self-disassembly or repair of the product outside of ROAMAN's special repair centers;
3. Lack of valid proof of purchase or the valid proof of purchase was altered;
4. Product damage due to voltage instability or voltage exceeding specified range;
5. Product exceeded warranty period or scope of warranty

After-sales service card:

Product model:

Date of purchase:

Customer name:

Contact information:

Fault description:

Warranty covers the unit and the charger (excluding the toothbrush head and the jet tip).

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