
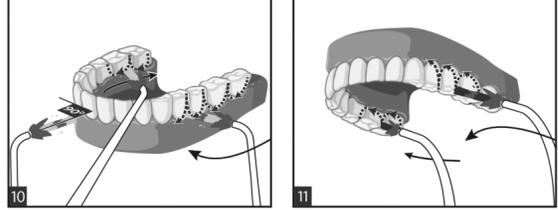
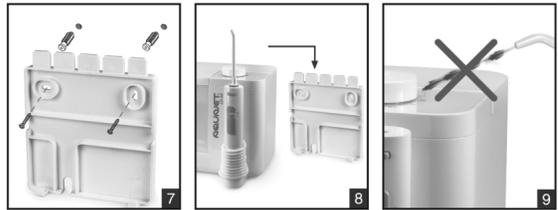
**ENG** Oral Irrigator AQUAJET LD-A9 Instruction Manual

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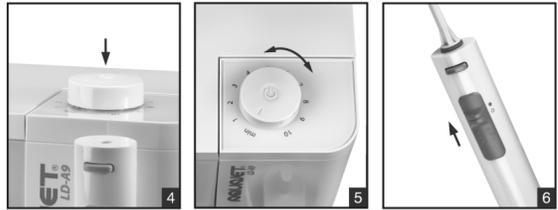**ENG**

**BEFORE USING THE AQUAJET LD-A9 ORAL IRRIGATOR (HEREINAFTER REFERRED TO AS THE IRRIGATOR, DEVICE, PRODUCT), PLEASE READ THIS INSTRUCTION MANUAL CAREFULLY AND RETAIN IT FOR FUTURE REFERENCE.**

**THE FOLLOWING SAFETY SIGNAL WORDS ARE USED THROUGHOUT THIS MANUAL:**

**⚠ WARNING!** INDICATES A POTENTIALLY HAZARDOUS SITUATION WHICH IF NOT AVOIDED, MAY RESULT IN INJURY OR PROPERTY DAMAGE.

**🚫 PROHIBITED!** INDICATES EXPRESSLY FORBIDDEN ACTIONS AND DEVICE HANDLING METHODS.

**📌 IMPORTANT!** INDICATES MANDATORY ACTIONS AND INSTRUCTIONS TO BE STRICTLY FOLLOWED TO ENSURE SAFE OPERATION.

Nº	Parts list	Nº	Parts list
1	Cover	5.3	LD-SA05 nozzle
2	ON/OFF button and water pressure regulation knob	6	Nozzle compartment
3	Water container	7	Water supply switch
4	Device body	8	Irrigator handle
5.1	LD-SA01 nozzle	9	Tube
5.2	LD-SA02 nozzle		

**GENERAL**

The oral mucosa is a complex organ of the human body that performs a wide variety of vital functions. The impaired blood circulation in the oral cavity, the presence of various dental structures, insufficient hygienic care lead to the emergence of various oral diseases, and, consequently, to a significant deterioration in the quality of life. Regular massaging effect of pulsating water jet on the soft tissues of the oral cavity stimulates blood flow, promotes disease prevention and therapy, which increases the natural protective properties of the oral mucosa to restore the protective, secretory, sensory, thermoregulatory and immune functions.

When it comes to dental care, most people tend to focus on "whiteness" of their teeth, neglecting the health of their gums. Meanwhile, failure to massage and clean the interdental space and gingival crevice (the area between the tooth neck and the gum line) leads to gradual tartar buildup, gum inflammation and persistent bad breath. As the inflammatory process develops, the gum tissue gradually moves away from the tooth, allowing bacteria to penetrate deeper, until it reaches the jawbone, whose structure is also subject to inflammatory destruction. Because of the bone tissue reduction, the gums also begin to recede. There is a visible lengthening of the teeth, and the smile no longer looks as aesthetically pleasing. The weakened support of the teeth by the gums and jawbones eventually leads to teeth loss. The space between your teeth is where bacteria accumulate and multiply. Bacteria live in the soft plaque and as a result of their activity acid is formed, which destroys the tooth enamel and leads to inflammation of the gum mucosa. Dental crowns and dentures are difficult to clean by regular methods, which causes pathogenic bacteria to accumulate between the edges of the crown and the gingival sulcus. This leads to the formation of tartar and inflammation of the oral mucosa.

**MEDICAL DEVICE PURPOSE**

The device is designed for massage and medicinal irrigation of the oral mucosa, which contributes to the prevention and therapy of periodontal disease, gingivitis, periodontitis.

**INDICATIONS FOR USE**

The irrigator forms pulsating water jet, which affects a particular area of the mouth, or the entire mouth, massaging the soft tissues. Massage of the mouth soft tissues restores and stimulates blood circulation, which in turn is essential for maintaining the health of the gums, teeth, and prolongs the life of the installed dental structures. Regular daily massage of the mouth soft tissues the most effective way to combat the natural loss of gums associated with age-related changes. In addition to the basic massage action, the device provides a thorough hygienic care of the entire oral cavity, including poorly accessible for a toothbrush/floss interdental spaces, gingival folds/pockets, dental structures such as dentures, bridges, crowns, braces. This prevents the development of pathogenic microflora and inflammatory processes in the contact areas of dentures and crowns with the oral mucosa, which is not achieved by other conventional means. Prevents formation of tartar, gingivitis and many periodontal diseases. Designed for use by adults and children.

- Contraindications:
- Contraindicated for patients with periodontal abscesses and oral ulcers.
  - Contraindicated for persons who cannot operate the device independently.

Side effects: No side effects have been identified.

**AQUAJET LD-A9 OPERATION BY CHILDREN**

**⚠ WARNING!** AQUAJET LD-A9 oral irrigator is an electrical device! This device can be operated by the child only under the supervision of an adult who reviewed this user manual.

The use of AQUAJET LD-A9 oral irrigator reduces the risk of dental diseases, facilitates oral care, and reduces the oral hygiene several times if the child wears braces. Also, the children often are not patient enough for a thorough tooth brushing, and the use of an irrigator helps maintain a superior oral hygiene.

The children can start using the oral irrigator since they begin using a toothbrush.

Despite the fact that the proper use of an oral irrigator cannot injure the mouth, always consult a pediatric dentist before the first use of an irrigator.

**CONTRAINDICATIONS**

**📌 Physician should review patients who previously require antibiotic premedication for dental treatment prior to introducing oral irrigation (bacteremia).**

Always follow the general safety precautions below when using electrical appliances, especially when children are present:

- ⚠ Always unplug this appliance immediately after using.**
- 🚫 Do not handle the plug with wet hands.**
- 🚫 Do not place or store product where it can fall or be pulled into a tub or sink, shower stall. Do not use while bathing.**
- 🚫 Do not place in or drop into water or other liquid.**
- 🚫 Do not reach for any appliance that has fallen into the water. Unplug it immediately.**
- 🚫 Never play a jet of water onto the gap between the device body and the ON/OFF button (Fig. 9). Otherwise, the liquid may get into the control unit and damage the device.**
- Before using, check the power cord for damage.**
- Do not use this appliance if the power cord or plug is damaged, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Return the appliance to an authorized service center nearby for examination and repair. Do not try to fix it yourself.**
- When the appliance is used by, on, or near children, close supervision is very necessary.**
- An appliance should never be left unattended when it is plugged in. Use the appliance only for its intended use as described in this manual. Do not use attachments or accessories not recommended by the manufacturer.**
- Keep the power cord away from hot or heated surfaces.**
- Remove the Water Container before filling. Fill with lukewarm water or appropriate oral rinse products only. Do not overfill.**
- Never drop or insert any object into any opening or hose.**
- Do not operate where oxygen or aerosol sprays are being used.**
- Do not use outdoors.**
- Do not modify this device without authorization of the manufacturer.**
- Do not replace the internal components until service has been performed by an authorized representative or manufacturer.**

**⚠ The device may only be used by children under the close supervision of adults.**

**All major maintains on the device must be performed by an authorized service center or distributor. The fuse is in the device, which should be changed by authorized service if damaged.**

**The power supply cord is replaceable only by professional maintenance personnel or manufacturer.**

**HOW TO USE**

**Placing Water Container.**

**🚫 Use of alcohol containing liquids, liquids with suspended matters, oily liquids, as well as other liquids not intended for irrigator application can entail irrigator failure!**

Open a cover and remove a water container. Fill the container with warm water or special solution.

**📌 When using for the first time, fill the container with warm water, press the On/Off button, set the jet pressure knob to the "min" position and slide the water supply switch up. Then move the knob gently several times to position "10" and back, so that water fills all internal cavity. The device is ready for first use when the water in the container is completely used up.**

Make sure that the temperature of the liquid in the container does not exceed 40°C. Do not pour hot solutions or water into the container, because the oral mucosa is very sensitive and can be damaged by hot solution. Put the container back in place (Fig. 1).

**Power Cord And Plug.**

The power cord is inserted into a special back-wall compartment (Fig. 2). Slide down the back cover and take out the power cord. Place the back cover to its original position, route out the power cord through a hole provided. Plug the power cord into a power outlet.

**Attaching/detaching nozzle to the irrigator handle.**

The device package includes nozzles with colored rims of different colors. Select an individually colored nozzle for each family member. Always use your individually colored nozzle only. This will allow full compliance with the requirements of personal hygiene.

Insert the nozzle into the irrigator handle until it clicks so that the nozzle guide aligns with the groove on the irrigator handle (Fig. 3a). Now, this device is ready for use.

To detach the nozzle, press the nozzle release button and pull the nozzle while holding it by the notches (Fig. 3b).

**Safe nozzle removal recommendations:**

1. If the nozzle does not detach well from the irrigator handle or slips in your wet hands, use a dry cloth. Wrap a light, dry cloth around the nozzle and detach it.
2. When removing the nozzle from the irrigator handle, DO NOT try turning the nozzle about its axis. Rotating the nozzle may damage it.

**📌 Additional nozzles can be purchased separately. Buy only AQUAJET LD-SA01, LD-SA02 and LD-SA05 nozzles for this device.**

**Power On Button.**

Press the power on button (Fig. 4). At that, green indicator will light up. For the first time use, please turn on the power on button for seconds until water coming out.

**Adjust Water Pressures.**

The jet pressure can be continuously and smoothly adjusted by turning the knob (Fig. 5). Initially, it is recommended to set the lowest pressure (the jet pressure knob is in min position), especially if your gums tend to bleed or if children are using the AQUAJET. Gradually increase the water pressure over time until finding the most comfortable setting for yourself.

Consult your professional dentists for recommendation of water pressure setting, particularly if the consumers have mouth disease.

**Start to Use.**

**⚠ Make sure that the temperature of the liquid in the container does not exceed 40°C.**

**The device's continuous operation should not exceed 5 minutes, this will be reminded by the light indicator, which will start blinking shortly before the end of this time. If you use the device longer or forget to turn it off, it will do so automatically (the device turns off after 5 minutes).**

Lean over the sink. Bring the nozzle tip to your teeth or gums. Start up water jet with the water switch located on the handle of irrigator (Fig. 6) by sliding it up. Keep your mouth slightly open to allow the water to run into the sink. Point the nozzle perpendicularly to the front teeth of the lower dental arch. In slow movements, for 4–5 seconds, clean the interdental space, holding the jet near the gums. Then, run the jet along the gum line at the base of the teeth to the next interdental space. Repeat these movements as for the first interdental space (Fig. 9), the dotted arrow indicates the place where the water jet hits, the solid arrow is the direction of nozzle movement. You can shut off the jet at any time by sliding the water supply switch on the irrigator handle. Do not direct the jet directly into the gum pocket, try to keep the jet perpendicular to the tooth line. Do the same procedure with all the interdental spaces on the outside and inside of the bottom teeth, one by one. Then repeat the procedure for the outside and inside of the upper dental arch (Fig. 10). If there are artificial prosthetic structures in the mouth, their places of contact with the gums, teeth or other dental appliances should be washed more carefully with the jet.

**Completing the Procedure**

Turn the device off by sliding the switch on the irrigator handle and turn off the power with the On/Off button on the main unit body.

**⚠ After use, always empty the Water Container, spray tube and nozzle entirely to prevent creating any bacteria from the standing water. Do not turn off the device until all water is out of the container.**

**WALL MOUNTING INSTALLATION**

Take off the rear panel by pushing it downward. Drill two bores in the wall and drive two plastic inserts into. Screw in the rear panel (Fig. 7). Take the power cord out of the storage box and fit it through the hole of rear panel. Do not press too hard to damage the power cord. Fix the device on the back cover (Fig. 8).

**GUARANTEE**

The purchased device is warranted for 2 years from the date of sale. Warranty obligations shall be formalized with a warranty card when the device is sold to a customer. Faults detected during the warranty period will be repaired free of charge by a specialized service center within 21 business days from the date of device delivery to the service center. The costs for the defective device transportation are borne by Little Doctor Europe Sp. z o.o. (door-to-door service), except for faults that are not covered by the warranty. For the addresses of the warranty service providers, see MANUFACTURER AND AUTHORIZED REPRESENTATIVES or go to [www.littledoctor.sg](http://www.littledoctor.sg).

The warranty does not cover the following: consumables and accessories for individual use (tubes, nebulizers, masks, mouthpieces, etc.); products with mechanical, thermal or chemical damage, or damage caused by violations of the rules specified in the operating instructions; products with signs of modification, opening and/or repair by an unauthorized service center (by an individual); products that have been connected to unrecommended power supplies or used with unrecommended consumables; products that have been operated with non-remedied defects; products that are defective through no fault of the manufacturer, such as: actions of third parties, natural phenomena and natural disasters, animals, insects, foreign objects or liquids entering the device, etc.; internal and external contamination, scratches, cracks, scrapes and other mechanical damage to the equipment occurring during its operation; products with illegible or altered serial numbers or a damaged warranty seal; batteries, glass parts, light bulbs, as well as operations such as adjusting, cleaning and other maintenance of the product.

The buyer has the right to replace the device with a serviceable one if an authorized service center finds a manufacturing defect that cannot be eliminated or after 3 unsuccessful repairs of the same malfunction. The party submitting the claim must properly prepare the device for shipment, protecting it from damage during transportation and enclosing a completed warranty card (with the seller's stamp, date of sale and seller's signature) and proof of purchase (invoice or receipt). The warranty does not cover damages resulting from inadequate protection of the device during transportation. The warranty period is extended for the time the device is under repair. Post-warranty repair of the device is performed at the user's expense.

**MAINTENANCE, STORAGE, REPAIR, AND DISPOSAL**

1. The device is to be sanitized after shipping and storage and after each use.

Disconnect the power cord from the power outlet before cleaning the device. Wipe the device outside surface, including the irrigator handle, with a cloth soaked in a solution with antimicrobial and detergent properties (soapy water, mild detergent, hydrogen peroxide, etc.), wringing it out beforehand. The nozzle compartment can be removed and the space under the nozzle compartment should also be cleaned periodically.

**⚠ Do not use rough cloths, brushes, or abrasives to clean the device.**

**Use only mild detergents. Never use solvents or alcohol, as they can cause malfunction or damage the device components.**

Wipe all treated surfaces dry with a soft cloth.

**📌 In the conditions of dental treatment and prevention facilities, nozzles shall not be disinfected to be passed on for use by other patients. After a patient treatment session is completed, the medical device is to be decontaminated while the accessories are to be disposed of.**

The container, nozzles, and nozzle compartment can be cleaned with mild detergents or in the dishwasher.

2. The device must be protected from direct sunlight and shocks.
3. Do not store or use the device close to heat sources and open flame.
4. Protect the device from contamination.
5. Do not allow the device to come into contact with corrosive solutions.
6. For any repairs please use authorized service centers only.
7. The month and year of production is indicated in the serial number after the letter A (where the first two digits indicate the year, the second two digits indicate the month). After the expected service life is over, periodically contact an authorized service center to check the device's condition and, if necessary, dispose of it per the local regulations. The manufacturer does not provide any disposal instructions.

**SPECIFICATIONS**

Model	LD-A9
Power source	~100-240 V / ~50/60 Hz
Power consumption	60 VA
Maximum water pressure	not less than 930 kPa
Water jet pulsation frequency	400–1800 pulses/min
Maximum continuous use time	5 min
Automatic shutdown	4 min 45 s ± 5 s

Water flow rate, min/max	≤250 ml/min / ≥330 ml/min
Complete set	Main unit, container, AQUAJET LD-SA01 nozzles (2 pcs.), AQUAJET LD-SA02 nozzle (1 pc.), AQUAJET LD-SA05 nozzle (1 pc.), panel for wall mounting, set of screws and wall plugs, operation manual, warranty card, consumer package (packaging).
Water Container Volume, no less than	500 ml
Operating conditions:	
Temperature	10°C to 35°C
Relative humidity	30-80% Rh
Atmospheric pressure	86-106 kPa
Storage and transportation conditions:	
Temperature	-20°C to +40°C
Relative humidity	15-95% Rh
Atmospheric pressure	50-106 kPa
External dimension	160 (H) x 170 (W) x 168 (D) mm
Net Weight	1200 ± 20 g
Year and month of manufacture	indicated in the bottom of the unit body in a serial number as AYYMMA9XXXX, where YY is the year, and MM is the month of manufacture

**Symbols explanation:**

- Type BF protection
- Class IP protection
- Manufacturer
- Serial number
- Importer
- Medical Device
- Compliance with the Directive MDR (EU) 2017/745
- Protection class IP
- For disposal, refer to your local regulations
- Authorized Representative in the EU
- Unique Device Identification
- Important: Read the instructions

The device manufacturing facility is ISO 13485 compliant. The device complies with the EU Directive MDR (EU) 2017/745.

The revision date of this Manual is indicated on the last page as XX-XXXX-YYMM-NN, where YY is the year, MM is the month, and NN is the revision number.

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

**TROUBLESHOOTING**

Problem	Possible reason	Solution
The device does not work	No electricity in a socket Poor contact between the plug and the socket	Try other sockets Unplug the power plug and plug it back in
Pump is functioning, but there is no water sprayed out	There is no water in the container Jet pressure knob is in the lowest position	Fill the container with water Slide the jet pressure knob upwards
Water pressure is too low	Incorrect position of the jet pressure knob Incorrectly installed water container	Adjust the pressure with the jet pressure knob Correctly install the water container
The water container leaks	The container rubber valve is loose	Rinse the valve under running warm water using a detergent
Adjusting the jet pressure changes the pulsation frequency and the compressor sound.	You can adjust the pressure by changing the pulsation frequency.	This is not a defect.
The compressor stops when the jet pressure is set to minimum.	If you turn off the water supply on the handle with the jet pressure set to minimum, the compressor may stop.	This is not a defect. The compressor will turn on when you start the water supply again.

If in spite of the above recommendations you cannot achieve the correct results, stop using the device and contact the maintenance organization. Do not try to repair the internal mechanism by yourself.

**IMPORTANT INFORMATION REGARDING ELECTRO MAGNETIC COMPATIBILITY (EMC)**

This Device manufactured by Little Doctor International (S) Pte. Ltd. confirms to IEC60601-1-2:2014 Electro Magnetic Compatibility (EMC) standard. Further documentation in accordance with this EMC standard is available at Little Doctor International (S) Pte. Ltd. at the address mentioned in this instruction manual or at [www.littledoctor.sg](http://www.littledoctor.sg). Refer to the EMC information for Device on the website.

**MANUFACTURER AND AUTHORIZED REPRESENTATIVES**

**Manufactured under control and for Little Doctor International (S) Pte. Ltd.** (7500A BEACH ROAD #11-313 THE PLAZA SINGAPORE 199591. Postal Address: Senggang Central PO Box 960 Singapore 915411).

**Manufacturer:** Little Doctor Electronic (Nantong) Co. Ltd. (No. 8, Tongxiang Road Economic & Technical Development Area, 226010 Nantong, Jiangsu, People's Republic of China).

Little Doctor Europe Sp. z o.o. (57G Zawila Street, 30-390 Kraków, Poland, phone: +48 12 2684746, 12 2684747, fax: +48 12 268 47 53, e-mail: [biuro@littledoctor.pl](mailto:biuro@littledoctor.pl)).

Little Doctor Europe Sp. z o.o., 57G Zawila Street Krakow 30-390 Poland.

For more information please visit [www.aquajet.sg](http://www.aquajet.sg)

**ADDITIONAL IRRIGATOR ACCESSORIES\***

- AQUAJET LD-SA01 Nozzle**
  - Made of plastic.
  - For individual use.
  - Number in the package – 2 pcs.
- AQUAJET LD-SA02 Nozzle**
  - Special subgingival nozzle with a flexible nozzle for cleaning and hygiene of difficult areas.
  - Made of plastic with a flexible nozzle
  - For individual use.
  - Number in the package – 2 pcs.
- AQUAJET LD-SA05 Nozzle**
  - Special orthodontic nozzle with a toothbrush for clearing dental brackets.
  - Made of plastic.
  - For individual use.
  - Number in the package – 2 pcs.

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