Model:
MB-102
MW-102
MI-102
MR-102

SMART+

Sonic toothbrush

USER MANUAL

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The logo for SMART+ features the word "SMART" in white, uppercase letters inside a teal rounded rectangle, followed by a white plus sign (+).

Sonic toothbrush

First of all, thank you for choosing our electric toothbrush. We hope you will be satisfied with the quality of our service. Please read this manual carefully to learn how to use the toothbrush properly.

Important: Check for any obvious damage to the toothbrush before each use. If there is any damage, please do not use the toothbrush. If you have had a mouth or gum operation in the last two months, please consult your dentist before using this toothbrush.

If your teeth bleed after using this toothbrush, please see your dentist.

The toothbrush is only for cleaning teeth and gums. Please don't use it for anything else. If you feel uncomfortable or in pain when using the toothbrush, stop immediately and ask for help from your doctor. This product is not suitable for children, infants or people who are unable to supervise themselves, unless they are supervised by someone in charge.

01 PRODUCT DESCRIPTION

Your dentist recommends replacing your power toothbrush head every 2-3 months.

02 DISPLAY SCREEN

Power On

Normal mode is the standard brushing mode, ideal for daily oral care. It provides thorough removal of plaque and debris, taking care of the health of teeth and gums. It's a versatile choice for any user.

Clean

The Whitening Mode in the sonic toothbrush is designed to effectively remove stains from the surface of teeth. Thanks to its intense vibration and appropriate rhythm, it provides a whitening effect with daily use. It is ideal for people who want a brighter smile.

White

Sensitive

Sensitive mode is designed for those with sensitive teeth and gums. It works with softer vibrations to provide gentle yet effective cleaning without irritation. Ideal for daily care of sensitive mouths.

Refresh

Refresh mode is an intense brushing mode, ideal for quick and effective oral refreshment. With powerful vibrations, it provides deep cleaning of teeth and removal of food debris and plaque in a short period of time. Recommended for those looking for an intense cleaning effect.

Splash Proof

The Splash Proof mode automatically activates the toothbrush when the tip is gently pressed against the teeth. When the tip is pushed away or the pressure stops, the toothbrush automatically turns off.

„OVERFORCE“ SYSTEM

The sonic toothbrush is equipped with the Overforce system, which protects teeth and gums from too much pressure. When the pressure is too high, a red LED light on the back of the toothbrush will flash. Simply reduce the pressure to ensure effective and safe brushing.

CHARGING THE TOOTHBRUSH

When the message “Low power” appears on the display, charge the toothbrush using the USB-C cable. Connect the toothbrush to the charger and let it charge until “Finish” appears on the display, indicating that charging is complete.


CHANGING THE LANGUAGE

1. Turn on the toothbrush.
2. Hold down the power button for 10 seconds until RESET appears on the display.
3. Briefly click the power button to select the language of the reset program.
4. After entering the reset program, select the language again by clicking the power button.
5. hold the power button for about 3 seconds until the display shows Done.

The toothbrush is now set to the selected language.

03 USING

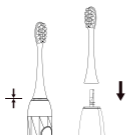
Step 1

Install the brush head firmly on the brush head shaft with a clearance of approximately 1mm between the brush head and the handle.


Step 2

Soak the toothbrush head and apply appropriate amount of toothpaste to the brush head.

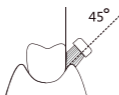
Step 3

Gently press the toothbrush head against the boundary between the teeth and gums.

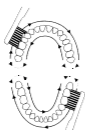Step 4

Press the button, select the proper brushing mode, and gently move the toothbrush head back and forth between the teeth to clean all parts of your teeth.

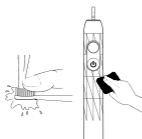**Step 5**

To prevent toothpaste from splashing, please turn off the toothbrush before taking it out of your mouth. After brushing, clean the toothbrush head with running water, shake off the water on the brush head, and wipe the handle clean with a soft cloth.

Rated voltage	3.7V
Rated power	2.5W
Waterproof level	IPX7
Charging time	8h
Battery capacity	2500mAh

04 THE WARRANTY BATTERY PROCESSING

CLEANING AND CARE

- Turn off the toothbrush, rinse the handle and cap under running water, then dry with a tissue.
- Store the toothbrush in a dry place.
- Replace the cap every 3 months or more often depending on usage.

WARRANTY DOES NOT COVER

- Damage caused by fortuitous events (e.g., flooding, fire).
- Mechanical damage due to the user's fault or by using unrecommended materials.
- Tampering with warranty seals or interference by persons outside the service center.
- Wear and tear resulting from natural use (e.g. decrease in battery capacity).
- Immersion in water for more than 30 minutes.

EQUIPMENT DISPOSAL

The symbol of the crossed-out trash garbage can indicates the prohibition of discarding the equipment with household waste.

Return the used equipment to an authorized recycling center to protect the environment and health.

For information on collection points, contact your local authorities or recycling centers.

The product complies with the
EU Directive on Electromagnetic
Compatibility and Electrical Safety

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
03 POUŽITÍ

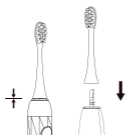
Step 1

Install the brush head firmly on the brush head shaft with a clearance of approximately 1 mm between the brush head and the handle.

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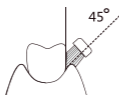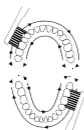
**Step 5**

To prevent toothpaste from splashing, please close the toothbrush before taking it out of your mouth.

After brushing, clean the toothbrush head with running water, shake off the water on the brush head, and wipe the handle clean with a soft cloth.

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